



# AIREGARD™ ES ENERGY SAVER HORIZONTAL LAMINAR FLOW CLEAN WORKSTATION

NU-240

## OPERATIONS & MAINTENANCE MANUAL



**OM0332**

Revision 3

March 2024



(115 VAC ONLY)

**Manufactured By:**

NuAire, Inc.

2100 Fernbrook Lane

Plymouth, MN 55447

Toll-Free: 1-800-328-3352

In Minnesota: (763)-553-1270

Fax: (763)-553-0459

# TABLE OF CONTENTS

|             |  |           |
|-------------|--|-----------|
| <b>1.0</b>  | <b>GENERAL DESCRIPTION</b>                               | <b>4</b>  |
| 1.1         | Safety Instructions                                      | 5         |
| 1.2         | Explanation of Symbols                                   | 5         |
| <b>2.0</b>  | <b>MODELS AND FEATURES</b>                               | <b>6</b>  |
| <b>3.0</b>  | <b>SHIPMENTS</b>   | <b>7</b>  |
| 3.1         | Damaged Shipments  | 7         |
| <b>4.0</b>  | <b>INSTALLATION INSTRUCTIONS</b>                         | <b>8</b>  |
| 4.1         | Location   | 8         |
| 4.2         | Set-Up Instructions                                      | 8         |
| 4.3         | Certification Testing Methods and Equipment              | 15        |
| <b>5.0</b>  | <b>OPERATING THE NU-240</b>                              | <b>16</b> |
| 5.1         | Aeromax™ Control System                                  | 16        |
| 5.2         | Operating Guidelines                                     | 21        |
| 5.3         | Ergonomics   | 21        |
| 5.4         | Cleaning Procedures                                      | 22        |
| <b>6.0</b>  | <b>GENERAL MAINTENANCE</b>                               | <b>23</b> |
| 6.1         | LED/Fluorescent Lamp Replacement                         | 23        |
| 6.2         | HEPA Filter Replacement                                  | 23        |
| 6.3         | Pre-Filter Replacement                                   | 23        |
| 6.4         | Airflow Control System Setup and Calibration             | 23        |
| 6.5         | HEPA Filter Leak Test                                    | 29        |
| 6.6         | Cleanliness Classification Test for Pharmacy Application | 29        |
| 6.7         | Main Control Board Description and Replacement           | 30        |
| <b>7.0</b>  | <b>OPTION PARAMETERS</b>                                 | <b>31</b> |
| 7.1         | Sync Function with Active Blower                         | 31        |
| 7.2         | Blower/Airflow Options                                   | 31        |
| <b>8.0</b>  | <b>ERROR MESSAGES AND TROUBLESHOOTING</b>                | <b>33</b> |
| <b>9.0</b>  | <b>POLYCARBONATE MATERIAL COMPATIBILITY</b>              | <b>34</b> |
| <b>10.0</b> | <b>REMOTE CONTACTS</b>                                   | <b>35</b> |
| 10.1        | Fan Relay  | 35        |
| 10.2        | Alarm Relay  | 35        |
| 10.3        | 15VDC Output   | 35        |
| <b>11.0</b> | <b>ELECTRICAL/ENVIRONMENTAL REQUIREMENTS</b>             | <b>36</b> |
| 11.1        | Electrical   | 36        |
| 11.2        | Operational Performance                                  | 36        |
| 11.3        | Light Exposure   | 36        |
| 11.4        | Installation Category: 2.0                               | 36        |
| 11.5        | Pollution Degree: 2.0                                    | 36        |
| 11.6        | Chemical Exposure  | 36        |
| 11.7        | EMC Performance  | 36        |
| <b>12.0</b> | <b>DISPOSAL AND RECYCLING</b>                            | <b>37</b> |

**Assembly Drawings**

ACD-06899.....Blower Brace Removal  
BCD-19015.....Base Stand Telescoping with Leg Levelers  
BCD-19016.....Base Stand Telescoping with Caster  
BCD-16521.....AEROMAX™ Front Panel  
BCD-19827.....Horizontal Airflow Filter Replacement  
BCD-20293.....Outlet option Instruction  
CD-000018.....Hydraulic Base Stand with Leg Levelers  
CD-000019.....Hydraulic Base Stand with Castors

**Electrical Schematics**

CD-000098.....Electrical Schematic 115V

**Purchase Specifications**

PSP-0018.....NU-240 Purchase Specifications

**AIREGARD™ ES Energy Saver Horizontal Laminar Flow  
Clean Workstation  
Model NU-240  
Manufactured By:  
NuAire, Inc.,  
Plymouth, Minnesota**

**1.0 GENERAL DESCRIPTION**











NuAire AIREGARD™ Horizontal Laminar Flow Clean Workstations utilize the newest technologies in laminar airflow design, materials, and manufacturing processes. The clean workstation can be used where clean airflow per ISO 14644-1 is required for the preparation of injectable drugs, IV solutions, tissue cultures, optics, microelectronics, etc. The clean bench **should not** be used for any work that involves biological agents assigned a level of Risk 1 through 4 as classified by the Centers for Disease Control (CDC), Atlanta, Georgia, since the horizontal flow offers no personnel protection against these agents. The clean bench should only be used to protect the product from contamination. The clean bench is optionally available with a base stand, placing the work surface at 30" (762mm) or 36" (914mm). A significant number of design innovations give the NuAire Laminar flow equipment superior performance qualities in airflow, lighting, noise levels and vibration.

## 1.1 Safety Instructions

These safety instructions describe the safety features of the AIREGARD™ Horizontal Laminar Flow Clean Workstation model NU-240. The workstation has been manufactured using the latest technological developments and has been thoroughly tested before delivery. It may, however, present potential hazards if it is not used according to the intended purpose or outside of operating parameters. Therefore, the following procedures must always be observed:

- The workstation must be operated only by trained and authorized personnel.
- For any operation of this unit, the operator must prepare clear and concise written instructions for operating and cleaning, utilizing applicable safety data sheets, plant hygiene guidelines, and technical regulations, in particular:
  - which decontamination measures are to be applied for the workstation and accessories.
  - which measures are to be taken in the case of an accident.
- Repairs to the device must be carried out only by trained and authorized personnel.
- Keep these operating instructions close to the unit so that safety instructions and important information are always accessible.
- Should you encounter problems that are not detailed adequately in the operating instructions, please contact your NuAire Representative or NuAire Technical Services.

## 1.2 Explanation of Symbols

| Symbol  | Explanation  |
|---|--|
|    | Safety alert symbol indicate a potentially hazardous situation which, if not avoided, could result in death or serious injury.                 |
|   | Safety alert symbol indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury.                 |
|  | CAUTION used without the safety alert symbol indicates a potentially hazardous situation which, if not avoided, may result in property damage. |
|  | Potential electrical hazard, only qualified person to access   |
|  | Biohazard  |
|  | Ground, Earth  |
|  | Lead Free  |
|  | Flammable Hazard   |
|  | Hazardous Gasses! Personal Protection Equipment Required.  |
|  | Chemical Hazard  |

## 2.0 Models and Features

NuAire's Model Number NU-240 designates the basic design series of AIREGARD™ Horizontal Laminar Flow Clean Workstation with the blower/motor located above the work surface (i.e. bench series). Model numbers are shown below.

### Model Number

**NU-240-330**

**NU-240-336**

**NU-240-430**

**NU-240-436**

**NU-240-530**

**NU-240-536**

**NU-240-630**

For detailed information reference [Purchase Specification PSP-0018](#).

## **3.0 Shipments**

NuAire takes every reasonable precaution to ensure that your AIREGARD™ Clean Workstation arrives without damage. Motor carriers are carefully selected and shipping cartons have been specially designed to protect your purchase. However, damage can occur in any shipment and the following outlines the steps you should take on receipt of a NuAire AIREGARD™ Clean Workstation to be sure that if damage has occurred, the proper claims and actions are taken immediately.


### **3.1 Damaged Shipments**

- 3.1.1** Terms are factory, unless stated otherwise. Therefore, it is important to check each shipment before acceptance.
- 3.1.2** If there is visible damage, the material can be accepted after the driver makes a notation on the consignee's copy of the freight bill. Then an inspection must be made to verify the claim against the carrier. This inspection is the basis of your filing the claim against the carrier.
- 3.1.3** If concealed damage is found, it is absolutely necessary to NOTIFY THE FREIGHT AGENT AT ONCE and request an inspection. Without this inspection, the transportation company may not accept a claim for loss or damage. If the carrier will not perform the inspection, an affidavit must be prepared stating that they were contacted on a certain date and that they failed to comply with the request. This, along with other papers in the customer's possession will support the claim.

## 4.0 Installation Instructions

### 4.1 Location

Within the laboratory, pharmacy, production area, etc., the ideal location for the clean workstation is away from personnel traffic lanes, air vents (in or out), doors and/or other sources of disruptive air currents. If drafts or other disruptive air currents exceed the face velocity of the HEPA filter, the potential exists for contaminated air to enter the work area of the workstation. Please note, however, that the work area has been constructed to minimize the effect of air currents by providing sidewalls that are an extension of the HEPA filter flow area.

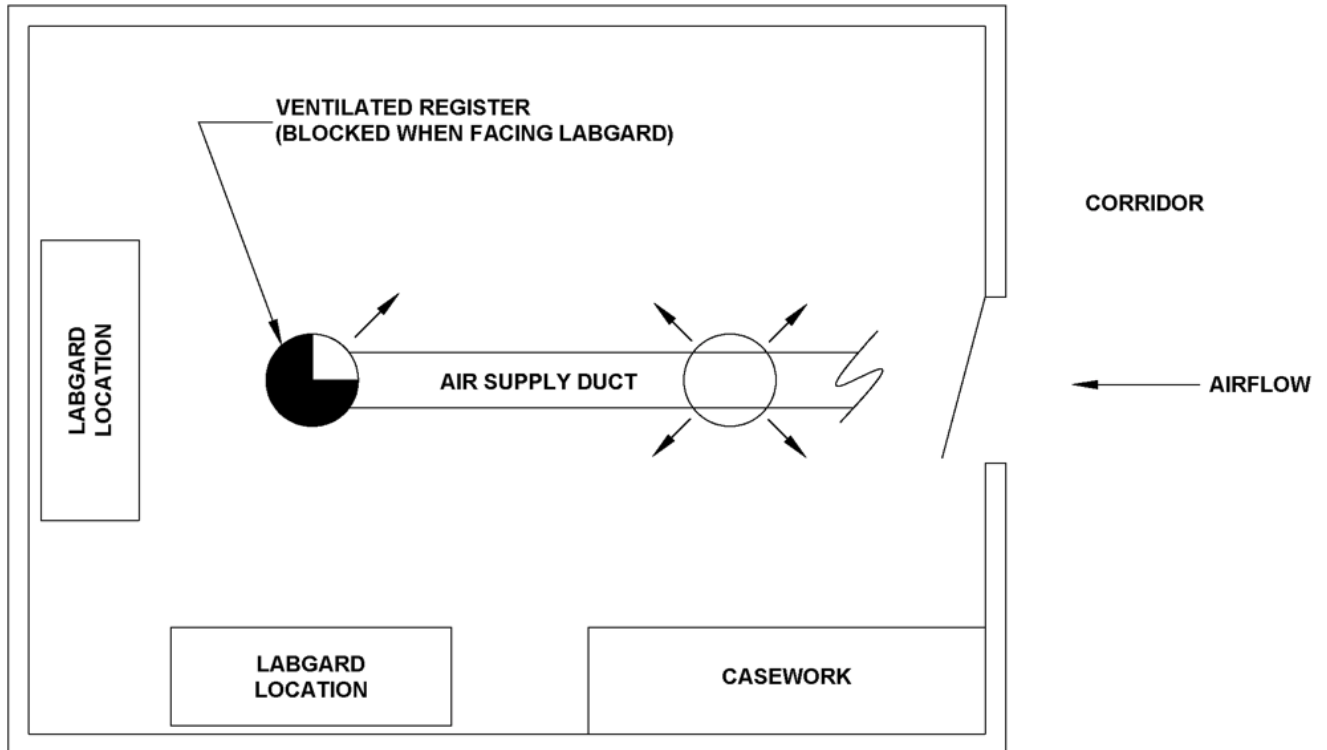

Where space permits, a clean six (6) inch (152mm) area should be permitted on each side of the workstation for maintenance purposes. The electrical outlet into which the workstation is connected should be readily accessible.

The workstation shall be positioned as not to obstruct the power supply outlet or the circuit breaker distribution panel.

### 4.2 Set-Up Instructions

Remove outer shipping protection (carton or crating). If the workstation is fastened to a base skid, it is usually the best procedure to leave the skid in place until the workstation is located in its approximate position to facilitate ease in handling. It can then be removed from the skid by removing the four bolts holding the workstation to the skid.

#### 4.2.1 Motor Shipping Bracket

The motor is securely fastened to the workstation during shipment to prevent damage to the blower housing mounting brackets. The procedure for removing the motor shipping bracket is shown in Drawing ACD-06899.

#### 4.2.2 Base Stand Assembly


The base stand is shipped K.D. in a separate carton and is assembled per the drawing accompanied with the unit. With a forklift or other suitable lifting device, lift the NU-240 between the top of the skid and the bottom of the workstation, slightly off of the floor and remove the four bolts holding the skid to the workstation (one in each corner). Now lift the workstation on top of the base and bolt the base stand to the workstation using four 4/20 x 3/4" (18mm) bolts and washers provided. Place the workstation in its desired location.



#### 4.2.3 Gas Service

NuAire doesn't recommend the use of natural gas within the clean bench, but if gas service is determined to be necessary for the application by the appropriate safety personnel, appropriate safety measures must take place.

Once the determination has been made by the appropriate safety personnel, the application of natural gas must be performed in accordance to national, state and local codes.

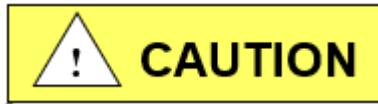

**IT IS ALSO STRONGLY RECOMMENDED THAT AN EMERGENCY GAS SHUTOFF VALVE BE PLACED JUST OUTSIDE THE CLEAN BENCH ON THE GAS SUPPLY LINE.**

All NuAire clean benches meet the safety requirements of UL and CSA for Laboratory Equipment. To comply with these safety requirements, NuAire uses only certified gas valves. In addition, if external piping is required, only black pipe is used for this application.

As previously stated NuAire doesn't recommend the use of natural gas within the clean bench and **ASSUMES NO RESPONSIBILITY FOR ITS USE. USE AT YOUR OWN RISK.** The Bunsen burner flame within the clean bench disrupts the laminar air stream, which must be maintained for maximum efficiency.



**IF THE PROCEDURE DEMANDS USE OF A FLAME, A BUNSEN BURNER WITH ON DEMAND IGNITION IS STRONGLY RECOMMENDED. DO NOT USE CONSTANT FLAME GAS BURNERS.**


During use, the Bunsen burner should be placed to the front of the workspace where resulting air turbulence will have a minimal effect.

#### 4.2.4 Plumbing Services

Ground key cocks with the type of service specified by the snap-in button on the handle, are located in the work zone. The ground key cocks are not recommended for pressures over 30 p.s.i. Reducing valves should be installed external to the workstation if necessary. Ground key cocks should never be used for oxygen service. A special needle valve for oxygen service is required and available upon request.

External connection is to 3/8 inch (10mm) FPT coupling through the sidewalls. Connection to plant utilities should be made with proper materials for the individual service and according to national and/or local codes. Observe all labels pertaining to the type of service and operating pressure.

|     |        |                              |            |      |      |
|-----|--------|------------------------------|------------|------|------|
| REV | CO     | DESCRIPTION                  | DATE       | DFTR | CHKD |
| B   | 000035 | UPDATED NOTES AND TITLEBLOCK | 05/30/2019 | LS   | BK   |

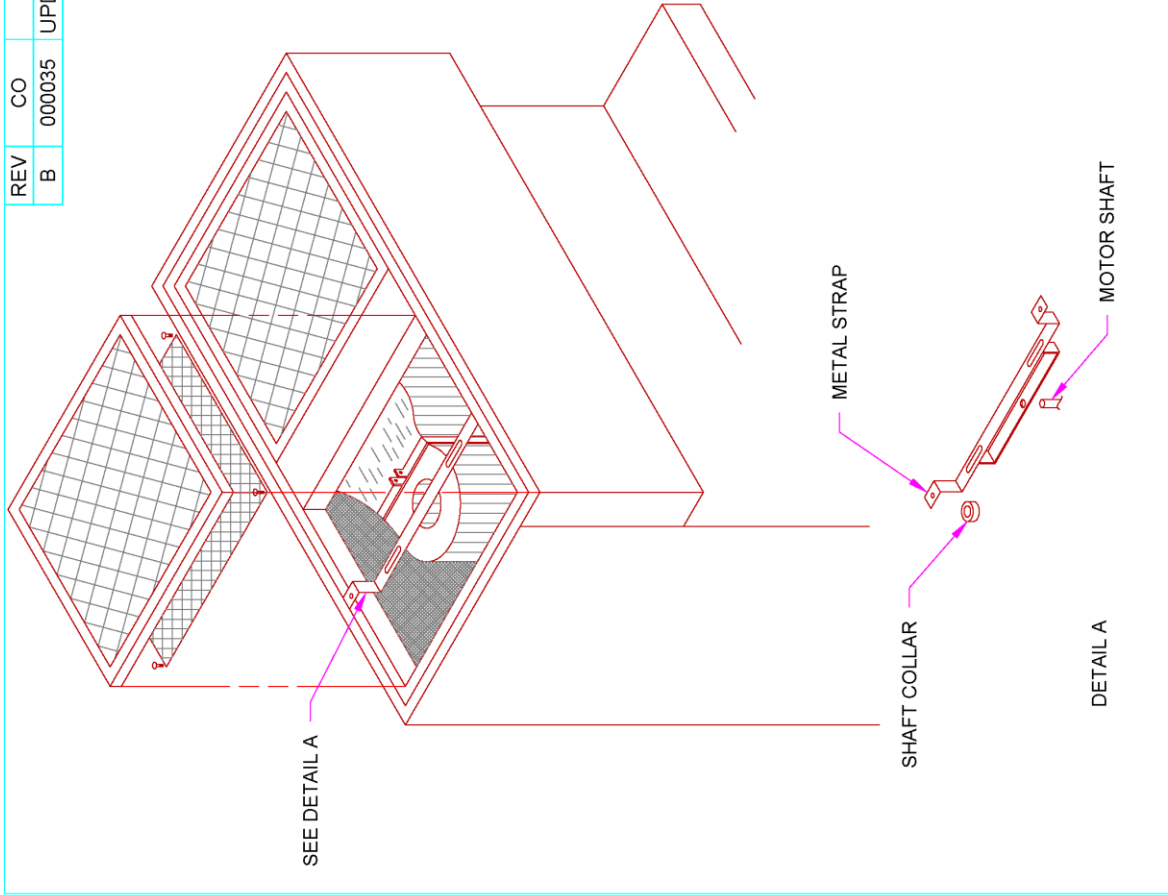

- STEPS:
1. REMOVE PRE-FILTER ON TOP OF UNIT ON BLOWER SIDE.
  2. REMOVE THE PRE-FILTER SCREEN SCREWS AND THE PRE-FILTER SCREEN.
  3. REMOVE THE THE 4 SCREWS ON THE METAL STRAP.
  4. LOOSEN THE SET SCREW ON THE SHAFT COLLAR.
  5. REMOVE THE STRAP AND COLLAR FROM THE UNIT.
  6. REPLACE THE PRE-FILTER SCREEN AND SCREW INTO PLACE.
  7. REPLACE THE PRE-FILTER.


NOTE: 1. THE METAL STRAP IS ONLY TO STABILIZE BLOWER SUSPENSION AND MOTOR AND MUST BE REMOVED FOR PROPER OPERATION.




|                |    |                      |         |
|----------------|----|----------------------|---------|
| TITLE          |    | BLOWER BRACE REMOVAL |         |
| DFTR           | PS | DATE                 | CHKD AS |
|                |    | 11/15/1988           |         |
| DRAWING NUMBER |    | ACD-06899            |         |
|                |    | SHEET                | 1 OF 1  |
|                |    | B                    |         |

| ITEM | QTY | PART NUMBER   | 3 FT. | 4 FT. | 5 FT. | 6 FT. | DESCRIPTION                          | REV | ECO   | DESCRIPTION      | DATE       | DFTM | CHKD |
|------|-----|---------------|-------|-------|-------|-------|--------------------------------------|-----|-------|------------------|------------|------|------|
| 1    | 2   | A-425-7373    |       |       |       |       | Lateral Biosafety Cabinet Latch      | C   | 13903 | UPDATED HARDWARE | 11/14/2018 | TH   | DSH  |
| 2    | 8   | A-500-395     |       |       |       |       | Basestand Leg Bushing                |     |       |                  |            |      |      |
| 3    | 1   | B-201-945-TAB | -01   | -02   |       |       | Upper Cross Member Panel Weldment 5" |     |       |                  |            |      |      |
| 4    | 2   | B-400-9562    |       |       |       |       | Basestand End Panel Weldment 28"     |     |       |                  |            |      |      |
| 5    | 4   | B-500-390     |       |       |       |       | Basestand Leg Extension              |     |       |                  |            |      |      |
| 6    | 1   | B-500-391-TAB | -01   | -02   | -03   | -04   | Cross Member Panel Weldment 5"       |     |       |                  |            |      |      |
| 7    | 4   | B-500-393     |       |       |       |       | Front Mounting Bracket               |     |       |                  |            |      |      |
| 8    | 4   | X-980443      |       |       |       |       | LEG LEVELER FOR 2" SQ TUBING         |     |       |                  |            |      |      |
| 9    | 8   | X-980493      |       |       |       |       | Acorn Nut 1/4-20 (Low Crown) SST     |     |       |                  |            |      |      |
| 10   | 8   | X-980494-03   |       |       |       |       | Nutsert 1/4-20                       |     |       |                  |            |      |      |
| 11   | 8   | X-980495      |       |       |       |       | Bolt HHCS 1/4-20 X 5/8 SST           |     |       |                  |            |      |      |
| 12   | 26  | X-980659-12   |       |       |       |       | Lock Washer 3/8 Internal Tooth SST   |     |       |                  |            |      |      |
| 13   | 8   | X-980781-07   |       |       |       |       | Washer 3/8 SAE Nickel                |     |       |                  |            |      |      |
| 14   | 16  | X-980783-07   |       |       |       |       | Washer 1/4 SST                       |     |       |                  |            |      |      |
| 15   | 34  | X-980794-02   |       |       |       |       | Lock Washer 1/4 Internal Tooth SST   |     |       |                  |            |      |      |
| 16   | 26  | X-980815-01   |       |       |       |       | Bolt 3/8-16 X 2-3/4 SST              |     |       |                  |            |      |      |
| 17   | 8   | X-980929-01   |       |       |       |       |                                      |     |       |                  |            |      |      |

Parts List

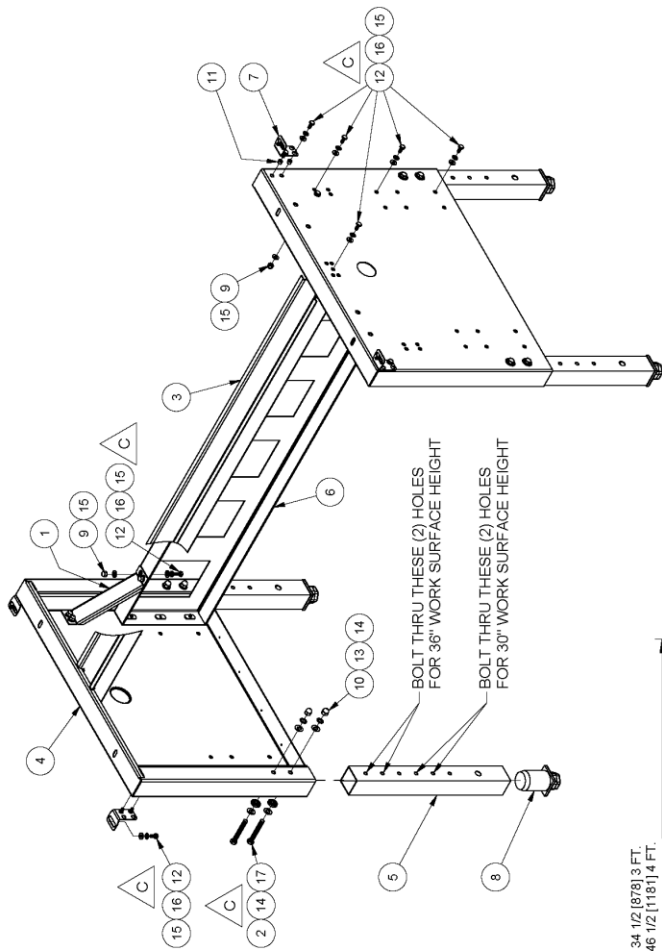

34 1/2 [878] 3 FT.  
 46 1/2 [1181] 4 FT.  
 58 1/2 [1487] 5 FT.  
 70 1/2 [1791] 6 FT.

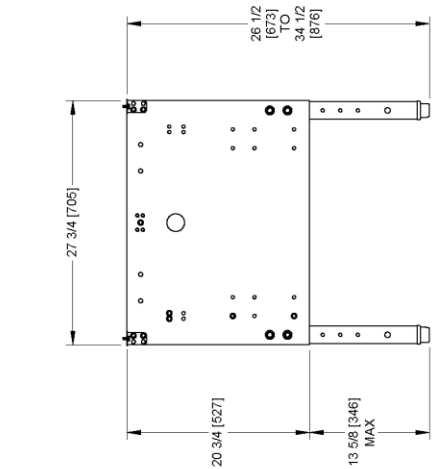


|            |
|------------|
| NU-201     |
| NU-200-133 |
| NU-200-134 |
| NU-200-135 |
| NU-200-136 |

|                       |           |           |       |   |    |   |
|-----------------------|-----------|-----------|-------|---|----|---|
|                       |           |           |       |   |    |   |
| TITLE ASSEMBLY 201    |           |           |       |   |    |   |
| BASESTAND TELESCOPING |           |           |       |   |    |   |
| DFTM                  | DATE      | CHKD      | SHEET | 1 | OF | 1 |
| TH                    | 4/10/2017 | MSS       |       |   |    |   |
| DRAWING NUMBER        |           | BCD-19015 |       |   |    |   |
|                       |           | C         |       |   |    |   |

SECTION A-A


| ITEM | QTY | PART NUMBER   | 3 FT. | 4 FT. | 5 FT. | 6 FT. | DESCRIPTION                               |
|------|-----|---------------|-------|-------|-------|-------|---|
| 1    | 2   | A-425-7373    |       |       |       |       | BRACE LATERAL, BSC                        |
| 2    | 8   | A-500-395     |       |       |       |       | BUSHING, BASE STAND LEG                   |
| 3    | 1   | B-201-945-TAB | -01   | -02   |       |       | W/LMT, UPPER CROSS MEMBER PANEL 5 FT.     |
| 4    | 2   | B-400-9562    |       |       |       |       | WELDMENT, END PANEL, BASESTAND (28")      |
| 5    | 1   | B-500-391-TAB | -01   | -02   | -03   | -04   | WELDMENT, CROSS MEMBER PANEL 5 FT.        |
| 6    | 4   | B-500-393     |       |       |       |       | FRONT MOUNTING BRACKET                    |
| 7    | 4   | B-500-507     |       |       |       |       | WELDMENT, LEG EXTENSION (CASTOR)          |
| 8    | 4   | B-520-2285    |       |       |       |       | BASESTAND, SAFETY ANGLE FOR 4 INCH CASTER |
| 9    | 4   | X-980040-08   |       |       |       |       | CASTOR, SWIVEL W/ BRACKET, 4", SST        |
| 10   | 16  | X-980490-06   |       |       |       |       | NUT HEX 5/16 - 18, FINISHED, SST          |
| 11   | 8   | X-980493      |       |       |       |       | ACORN NUT, 1/4-20                         |
| 12   | 8   | X-980494-03   |       |       |       |       | NUT, ACORN 3/8-16 18-8 SST                |
| 13   | 8   | X-980495      |       |       |       |       | 1/4-20 NUTSERT                            |
| 14   | 26  | X-980659-12   |       |       |       |       | SCREW, HHCS, 1/4-20 X 5/8 18-8 SST        |
| 15   | 16  | X-980667-01   |       |       |       |       | 5/16-18 X 3/4 BOLT                        |
| 16   | 8   | X-980781-07   |       |       |       |       | 3/8 LOCK WASHER                           |
| 17   | 16  | X-980783-07   |       |       |       |       | WASHER FL 3/8 SAE NKL                     |
| 18   | 32  | X-980790-06   |       |       |       |       | WASHER FL 5/16 SAE SST                    |
| 19   | 34  | X-980794-02   |       |       |       |       | 1/4 SST FLAT WASHER                       |
| 20   | 16  | X-980795-05   |       |       |       |       | WASHER LK 5/16" SPLIT, 18-8 SST           |
| 21   | 26  | X-980815-01   |       |       |       |       | 1/4 INTERNAL TOOTH LOCK WASHER            |
| 22   | 8   | X-980929-01   |       |       |       |       | BOLT 3/8-16 X 2-3/4 SST                   |

Parts List

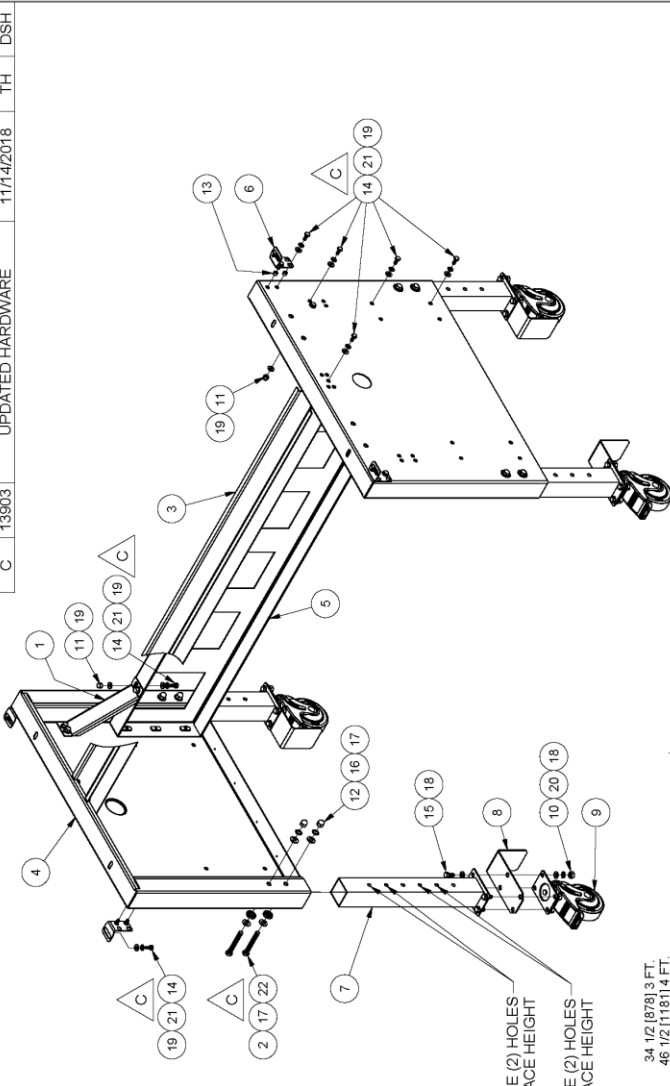

BOLT THRU THESE (2) HOLES  
FOR 36" WORK SURFACE HEIGHT

BOLT THRU THESE (2) HOLES  
FOR 30" WORK SURFACE HEIGHT

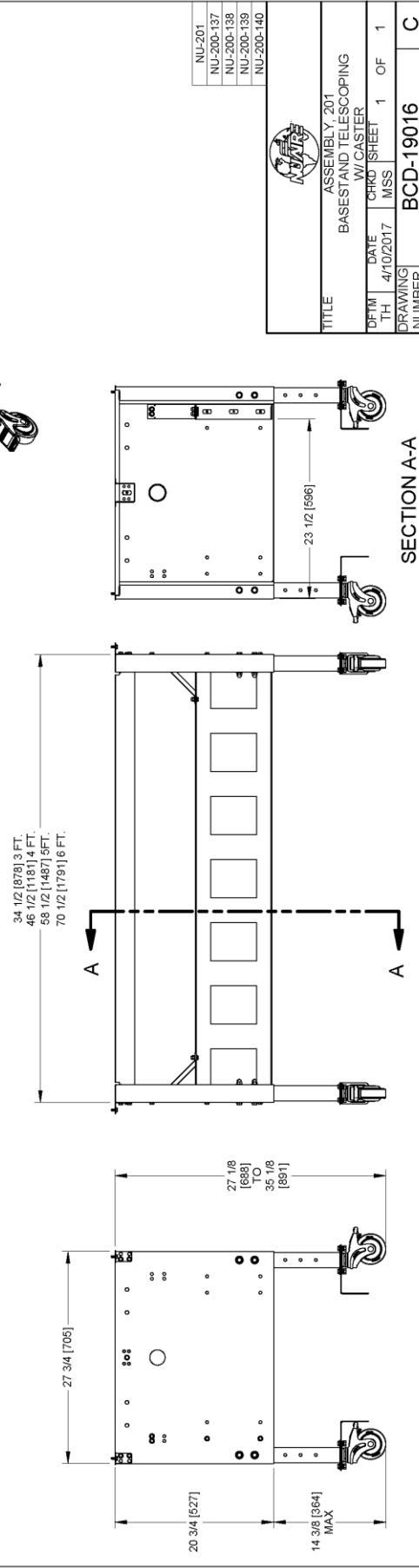


SECTION A-A

|     |       |                  |            |      |      |
|-----|-------|------------------|------------|------|------|
| REV | ECO   | DESCRIPTION      | DATE       | DFTM | CHKD |
| C   | 13903 | UPDATED HARDWARE | 11/14/2018 | TH   | DSH  |



  

|                |   |
|----------------|---|
|                |   |
| TITLE          | ASSEMBLY 201<br>BASESTAND TELESCOPING<br>W/CASTER |
| DFTM           | DATE  |
| TH             | 4/10/2017   |
| CHG            | SHEET 1 OF 1                                      |
| DRAWING NUMBER | BCD-19016   |
|                | C   |

| REV | CO     | DESCRIPTION            | DATE       | DFTM | CHKD |
|-----|--------|------------------------|------------|------|------|
| A   | 000132 | RELEASED TO PRODUCTION | 10/29/2019 | TH   | BK   |

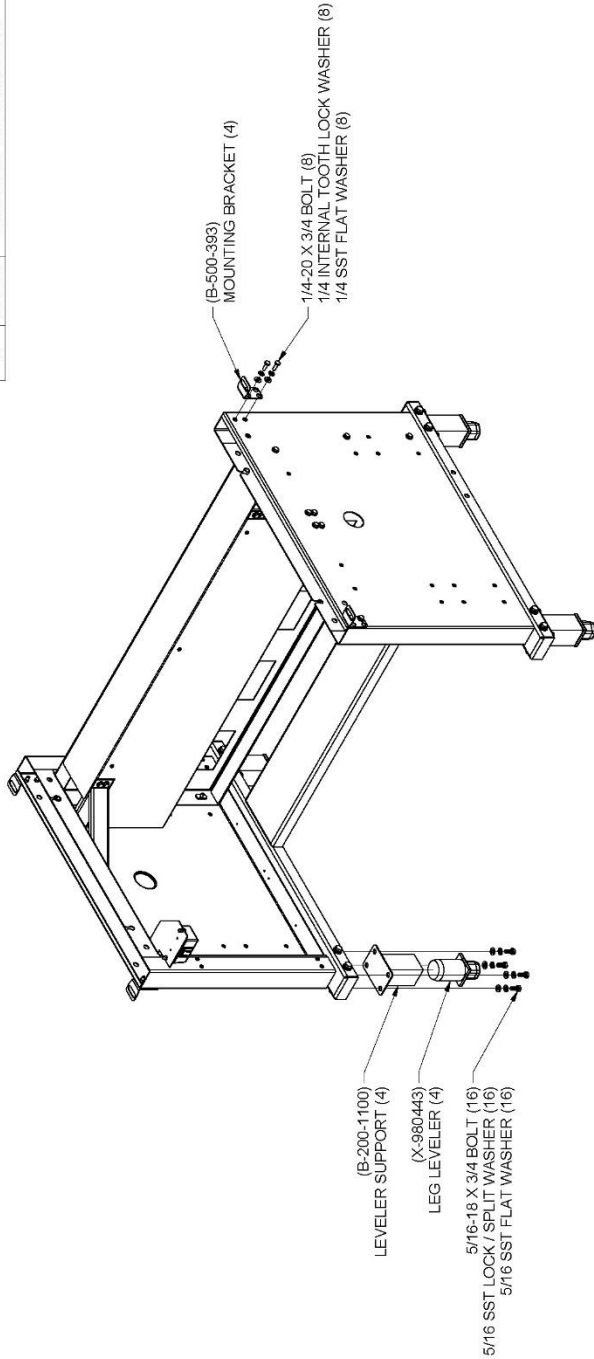

34 1/2 [876] 3 FT.  
 46 1/2 [1181] 4 FT.  
 58 1/2 [1486] 5 FT.  
 70 1/2 [1791] 6 FT.

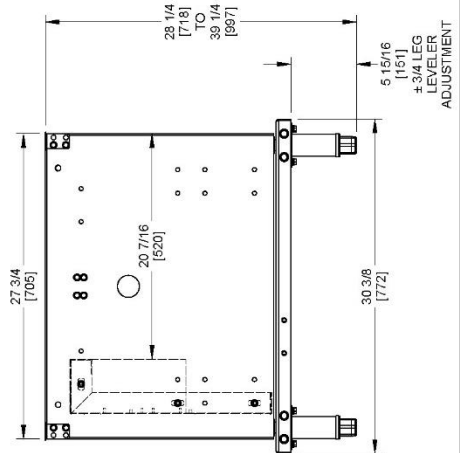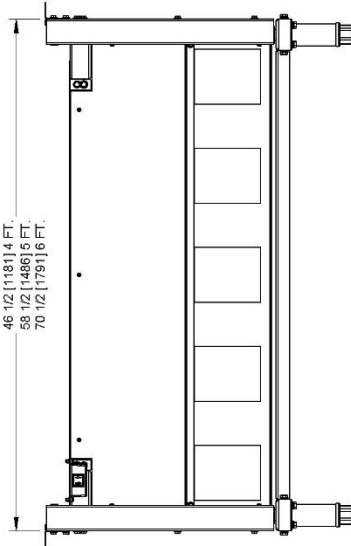



|            |
|------------|
| NU-200-147 |
| NU-200-148 |
| NU-200-151 |
| NU-200-152 |
| NU-200-155 |
| NU-200-156 |
| NU-200-159 |
| NU-200-160 |



TITLE  
 140/240 HYDRAULIC BASE  
 STAND ASSEMBLY W/LL

| DFTM           | DATE       | CHR | SHEET     | 1 | OF | 1 |
|----------------|------------|-----|-----------|---|----|---|
| TH             | 10/29/2019 | BK  |           |   |    |   |
| DRAWING NUMBER |            |     | CD-000018 |   |    |   |
|                |            |     | A         |   |    |   |

| REV | CO     | DESCRIPTION            | DATE       | DFTM | CHKD |
|-----|--------|------------------------|------------|------|------|
| A   | 000132 | RELEASED TO PRODUCTION | 10/29/2019 | TH   | BK   |

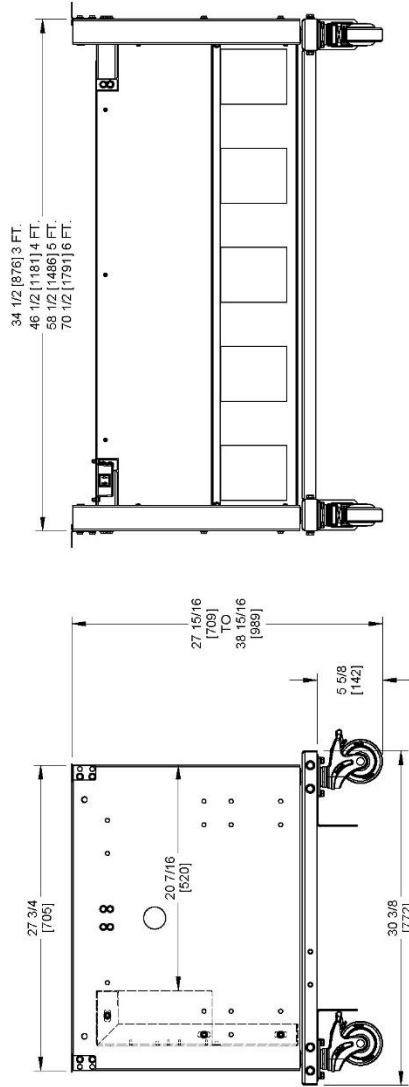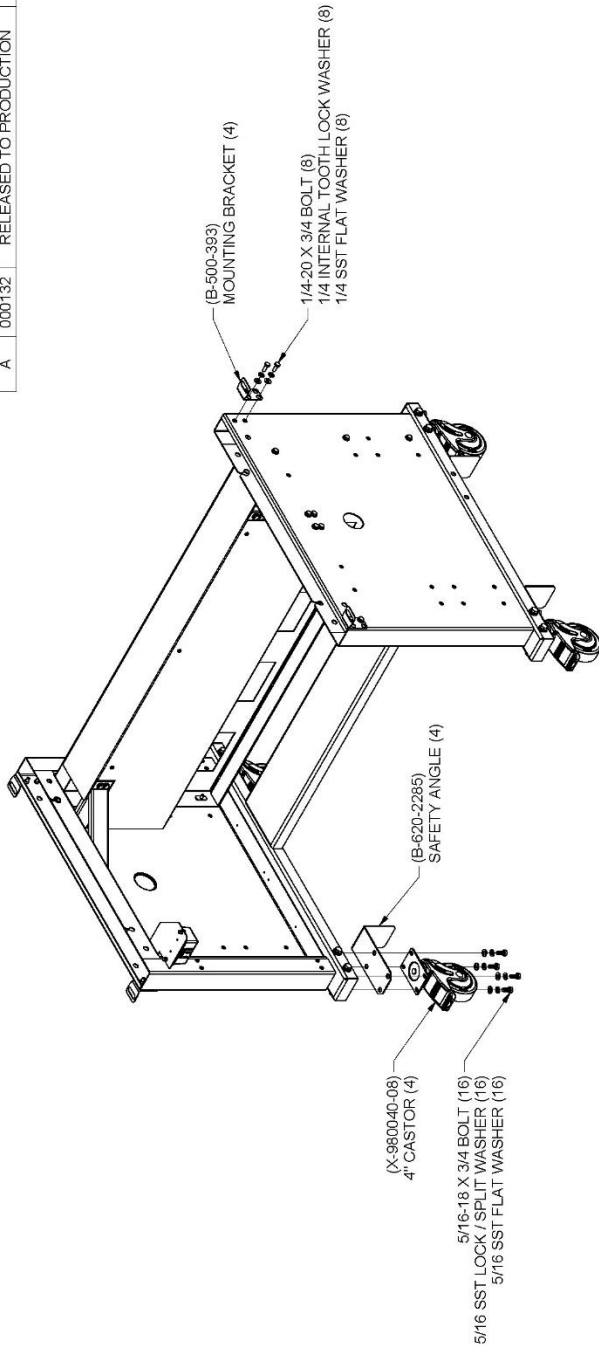

|            |
|------------|
| NU-200-149 |
| NU-200-150 |
| NU-200-153 |
| NU-200-154 |
| NU-200-157 |
| NU-200-158 |
| NU-200-161 |
| NU-200-162 |



140/240 HYDRAULIC BASE  
 STAND ASSEMBLY W/  
 CASTORS


| DFTM | DATE       | CHKD | SHEET | OF |
|------|------------|------|-------|----|
| TH   | 10/29/2019 | BK   | 1     | 1  |

DRAWING NUMBER: **CD-000019**

A

#### 4.2.5 Electrical Services

The Horizontal Laminar Flow Clean Workstation may be "hardwired" (optional) or connected via an electrical power cord, which is standard. The unit requires 115 VAC, 60 Hz, single phase. (Current rating varies per workstation size, reference Electrical/Environmental Requirements).

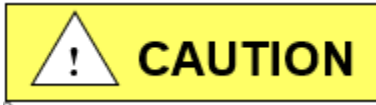

It is recommended that power to each work station (power is defined by each work station's power cord), whether individual or joined work stations, be on its own branch circuit, protected with a circuit breaker at the distribution panel.

**NOTE: THIS UNIT MAY CONTAIN ELECTRONIC BALLASTS FOR FLUORESCENT LIGHTING. ELECTRONIC BALLASTS OPERATE WITH HIGH INRUSH CURRENT. IT IS NOT RECOMMENDED TO USE THIS PRODUCT WITH GROUND FAULT CIRCUIT INTERRUPTERS (GFCI'S) BECAUSE THE BALLASTS MAY CAUSE THE GFCI TO TRIP.**

If a "hardwired" (conduit) connection is desired, an electrical junction box is provided on the top of the workstation with a removable cover. All wiring connections to the junction box should be done according to the National Electrical Code and/or local codes by a qualified electrician.

#### 4.2.6 Final Assembly

The exterior surfaces and viewing glass are easily cleaned with any mild household detergent cleaner using a soft cloth. Harsh chemicals, solvent-type cleaners and abrasive cleaners should not be used. See section 5.4 for cleaning procedures. Turn the workstation on and let it operate for 5 minutes before using it as a clean bench.

### 4.3 Certification Testing Methods and Equipment

After installation and prior to use, NuAire recommends that the workstation be recertified to factory standards. At a minimum, the following tests should be performed.

1. HEPA filter media
2. Filter frame leak test
3. Airflow velocities

The testing methods and equipment required are specified on the factory inspection report included with this manual.

**NOTE: IT IS RECOMMENDED THAT THESE TESTS BE PERFORMED BY A QUALIFIED TECHNICIAN WHO IS FAMILIAR WITH THE METHODS AND PROCEDURES FOR CERTIFYING CLEAN BENCHES.**

**NOTE: AFTER THE INITIAL CERTIFICATION, NUAIRE RECOMMENDS THAT THE WORKSTATION BE RECERTIFIED AT A MINIMUM ON AN ANNUAL BASIS AND AFTER EVERY FILTER CHANGE OR MAINTENANCE ACTION OR ANY TIME THE OPERATOR FEELS IT IS NECESSARY.**

**NOTE: SCANNING THE HEPA FILTER SEAL CANNOT BE DONE BY ONLY REMOVING THE DIFFUSER SCREEN. HEPA FILTERED AIR MUST BE DIRECTED INTO THE DEAD-AIR SPACE TO FLUSH THE AREA WITH PARTICLE FREE AIR WHILE SCANNING FOR LEAKS. THE FRAME ITSELF IS UNDER NEGATIVE PRESSURE AND SCAVENGER SLOTS, IF PRESENT, HELP FLUSH OUT THE AREA WHEN THE DIFFUSER IS INSTALLED.**

Note that NuAire Clean Workstations provide premium performance and superior reliability. However, protection of personnel and products is of utmost importance, and it is highly recommended that certifications are completed as stated to ensure consistent conformance to factory standards.

## 5.0 Operating the NU-240

### 5.1 Aeromax™ Control System

#### 5.1.1 Overview

The Aeromax™ control system is designed to service the control requirements of the AireGard™ ES NU-240. The Aeromax™ control system consists of an electronic module that will perform the following functions:

- Easy user interface via LED's and function keys
- Control blower via solid state switch.
- Control lights via solid state switch.
- Disable audible alarm switch with ring back function.
- Control blower DC ECM motor/blower with solid-state DC Motor Controller that provides automatic compensation for line voltage variances.
- Monitor and display airflow system performance via PresurFlow™ monitor.

The AireGard™ ES NU-240 offers the latest digital microprocessor design technology for improved performance and safety. The Aeromax™ control system integrates a digital pressure sensor (PresurFlow™) to monitor the airflow performance. The Aeromax™ control system also integrates a DC ECM motor/blower controller that provides automatic compensation for both filter loading and line voltage variances. There is additional on/off control of blower and light. All the above functions are shown in a system block diagram (see figure 1).

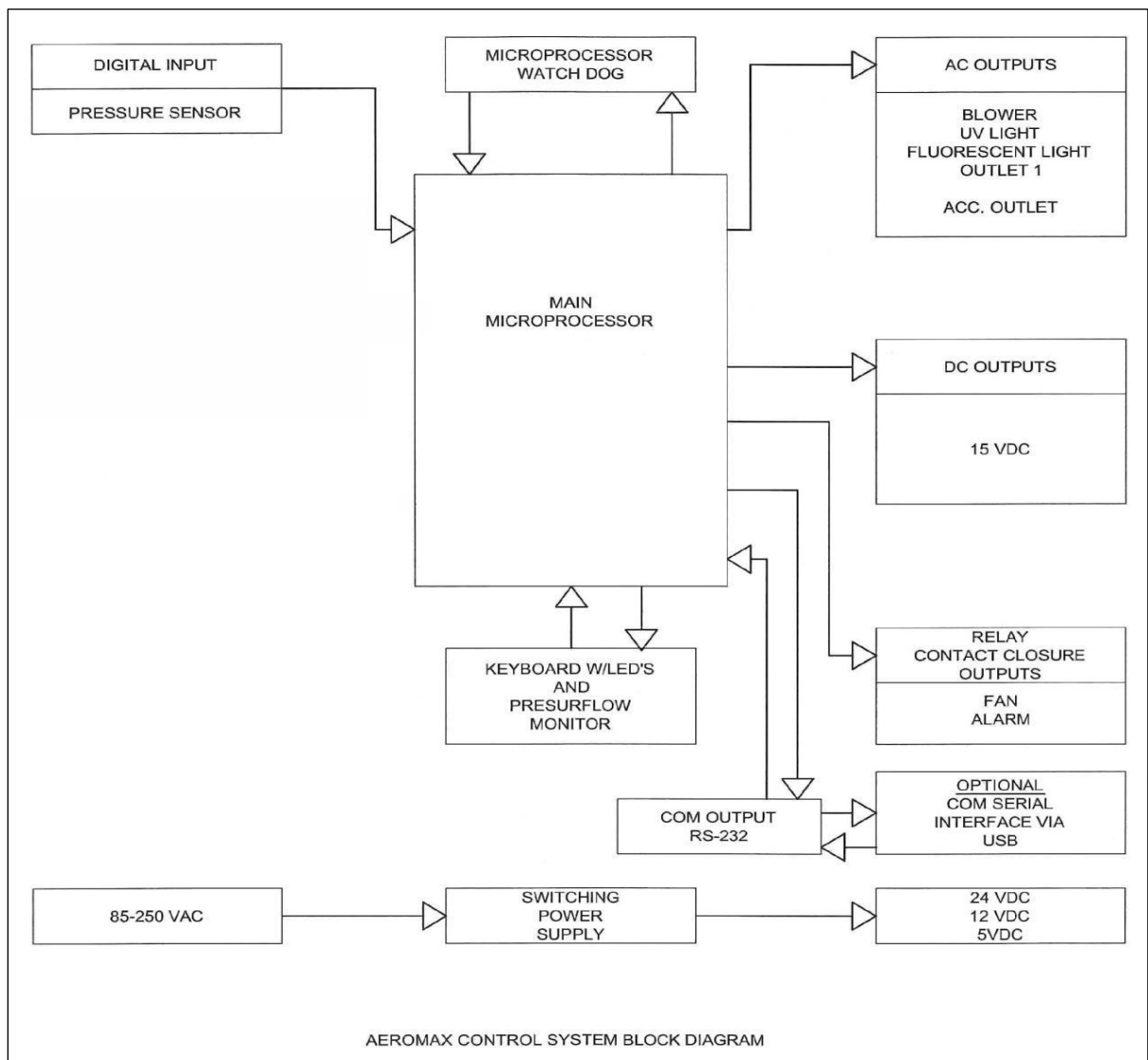


Figure 1



| REV | ECO   | DESCRIPTION   | DATE      | DFTM | CHKD |
|-----|-------|---------------|-----------|------|------|
| C   | 13471 | UPDATED TITLE | 8/25/2017 | TH   | BP   |



|                     |            |      |        |
|---------------------|------------|------|--------|
|                     |            |      |        |
| AEROMAX FRONT PANEL |            |      |        |
| DFTM                | DATE       | ORIG | SHEET  |
| U.S.                | 12/16/2013 | BP   | 1 OF 1 |
| DRAWING NUMBER      |            |      | C      |
| BCD-16521           |            |      |        |

### 5.1.2 Front Panel

The control system front panel contains the following functions described in detail (see Drawing BCD-16521).

- **Blower Keys**  
The blower key controls the ON/OFF power to the blower.  
LED above key indicates:
  - Full green for blower on
  - Blinking green for blower pending
  - Full red for blower alarm
- **Hidden Key**  
The hidden key is located just above the blower LED indicator centered in the airflow symbol. The hidden key is used for various functions including the blower password 3 key sequence if the option is activated.
- **LED Light Key**  
The LED light key controls the on/off power to the LED light.  
LED above the key indicates full blue for LED light on.
- **Ultraviolet (UV) Light Key (Diagnostic Key Use ONLY. UV Light feature not available on NU-240 models.)**  
The UV light key controls the on/off power to the UV light (Diagnostic Key Use ONLY).  
LED above indicates full yellow for UV light on.
- **Outlet Keys (Special Option Feature ONLY)**  
The outlet key controls the ON/OFF power to the outlets (Special Option Feature ONLY).  
LED above indicates full blue for outlets on.
- **Red Alarm LED**  
The red alarm LED will indicate any alarm condition and remain indicating until the alarm condition is cleared.  
**Audible Alarm Silence**  
The audible alarm silence key allows user interaction to silence an audible alarm for a period of 15 minutes. After 15 minutes if the alarm condition still exists, the audible alarm will again sound. The audible alarm silence key also is used to exit all Aeromax™ user interaction menus.
- **Arrow Adjustment Keys**  
The arrow adjustment keys allow user interaction for various functions.

### 5.1.3 Aeromax™ Control System Power

After the AireGard™ ES NU-240 is plugged into the appropriate facility line power the control system will power up. The control panel will also indicate the power up status by blinking the red alarm LED. Pressing any key will acknowledge the power up status and turn off the blinking red alarm LED.

If a power interruption occurs, all control system functions, calibrations and parameters will be maintained and continue upon restoration of power. Just as the initial power up, the red alarm LED will blink to indicate power up status.

### 5.1.4 Standby Mode

When the NU-240 is not in use any of the function keys except the blower that initiates run mode may be turned on and off in standby mode.

### 5.1.5 Run Mode

Any time the blower run key is pressed with the window at its correct operational height, the Run Mode screen will be initiated. While in Run Mode one LED segment will be illuminated to indicate cabinet airflow status. Wait for at least one minute after starting the blower for the airflow to stabilize before starting any work.

Once stabilized, the single illuminated LED will indicate the following:

- Green: Unit airflow is within safe working parameters.
- Yellow: Unit airflow is within safe working parameters, but near unsafe limits. Consider recalibration/recertification.
- Red: Unit airflows are not safe and unit should not be used. Check airflow values and recalibrate/recertify cabinet. See section 7.1.

**Note: When first setting the motor control function parameter (section 6.4.2) the PresurFlow™ system will enter an approximate 3 minute warm up period. If the blower is turned on during this period the PresurFlow™ LED indicators will blink and indicate the following sequence:**

- 1<sup>st</sup> minute – Left and right Red LED's will blink.
- 2<sup>nd</sup> minute – Left and right Green LED's will blink.
- 3<sup>rd</sup> minute – Center 3 Green LED's will blink.

Once the warm-up period is complete, only one LED will indicate cabinet airflow status.

### 5.1.6 Standby/Run Mode Alarms

If present, standby/run mode alarms will be both visual and audible, the red alarm LED oval will turn on. Audible alarms will produce an alarm tone for 30 seconds, then ring back for 2 seconds of every 5 seconds. Pressing the alarm silence key will silence the audible alarm for 15 minutes initially then will start the ring back function again.

The list below represents alarm types and their respective priority from the highest to lowest priority.



- 1) New Firmware Loaded
- 2) Internal Board Failure
- 3) Power on Reset
- 4) Airflow Pressure Alarm
- 5) Blower RPM Failure


**Note: The above messages are described in greater detail in section 7.**

### 5.1.7 Operator Accessible Functions

#### 5.1.7.1 Access and Navigation

To access the operator accessible functions,


- Press and hold the  key, then enter the 3 key sequence for the desired function, then release the  key and follow each instruction set.

**Note: Pressing the  key at any time will abort and exit the process without saving any changes made. Pressing the hidden key will accept all changes and exit.**


#### 5.1.7.2 Auto Timer Duration

Auto timer duration timers are countdown timers for the functions displayed once time is entered into a function. The timer will begin to countdown upon the start of that function (i.e. press light key to start timing the light). The LED indicator above the function key will start to blink indicating the timer function. If the LED indicator is full on, no timer function is present. As the timer expires the function will turn off.


- Select auto timer duration function
  - Outlets (if equipped)

Press and hold  key, then press hidden – outlet – outlet keys sequentially.

LED indicator above outlet will blink fast. Adjust desired time as described below.
  - Lights


Press and hold  key, then press hidden – light – light keys sequentially.

LED indicator above light will blink fast. Adjust desired time as described below.
  - UV Light (Diagnostic key use only. UV Light feature not available on NU-240 models.)

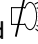Press and Hold  key, then press hidden – UV light – UV light keys sequentially.

LED indicator above UV light will blink fast.

Adjust desired time as described below.
- Low Flow Blower (Diagnostic key use only. Function not available on NU-240 models.)

Press and hold  key, then press hidden-blower-blower keys sequentially.

LED indicator above Blower key will blink fast. Adjust desired time as described below.

  - Adjust countdown time.
  - Press ↑ or ↓ keys to adjust time.

- Time will change in 15-minute increments as shown on the PresurFlow™ LED segments below.

| Time in Increments                | LED Display |
|-----------------------------------|-------------|
| 15 minutes for Red LED            |             |
| 30 Minutes for Yellow LED         |             |
| 1 hour for green LED              |             |
| Represents 8 hours (maximum time) |             |

- Press hidden key to accept time and exit.
- Press key at any time to abort and exit.

### 5.1.7.3 Blower Password

The blower on/off password allows the cabinet user to place a 3 key sequence requirement to turn the blower on or off.

The 3 key sequence for the blower password will be a combination of the hidden and blower keys.

- Select blower password  
Press and hold key, then press hidden – blower – hidden keys sequentially. Red LED indicator above blower will blink fast.
- Select password  
Press ↑ or ↓ key to scroll through the code choices below,

|  |                          |
|--|--------------------------|
|  | No Password              |
|  | B-B-H password           |
|  | B-H-B password (default) |
|  | B-H-H password only      |
|  | H-B-B password           |
|  | H-B-H password           |
|  | H-H-B password           |

- Press hidden key to accept time and exit.
- Press key at any time to abort and exit.

**Note:** If the required blower password option is selected in the blower airflow option menu (see section 7.2.2). Then the “No password” choice above is not available and the default remains B-H-B.

## 5.2 Operating Guidelines

Operate the laminar flow workstation continuously. The unit will then remain in its initially clean condition. If, for any reason, the unit is turned off, turn the unit on and permit to operate for 5 minutes before resuming operations.

Allow only essential items in the work station. Objects should not be placed between the HEPA filter and any point where the clean environment must be maintained. New items introduced into the work area should be placed downstream of items already in the work zone for several minutes to allow contaminants to flush off. Note that plastic parts may carry a static charge which may require special handling in order to remove contaminants.

Particular care must be exercised in placing equipment within the work space. Where possible, equipment should be placed on perforated platforms to allow air movement under as well as around the object.

All work should be performed with the operator's hand or head downstream of the critical process points. Unnecessary movement with the work station should be kept to a minimum.

If the workstation is used in a manner not specified by NuAire, the protection provided by the equipment may be impaired.

### 5.2.1 Operating Sequence

#### A. Start Up

Turn on workstation blower and lights, check air intake ports of the workstation to make sure they are unobstructed. The workstations are provided with gauges which indicate pressure differentials across the filters. They indicate when to replace the filters, dependent upon the blower fan capacity. Blower speed must only be readjusted by qualified maintenance technicians.

#### B. Good procedure includes the decontamination or wipe down of workstation surfaces with chemical disinfectant before work commences.

#### C. Allow blowers to operate for a minimum of 5 minutes before aseptic manipulations are begun in the workstation. An additional advantage is obtained from purification (filtration) of the room air circulated through the equipment. Because of the characteristic contributed to the quality of the laboratory environment, some owners leave them in operation beyond the time of actual use.

#### D. Minimize Room Activity - Activity in the room itself should be held to a minimum. Unnecessary activity may create disruptive air currents, as well as interfere with the work of the operator. A person walking past the front of the workstation can cause draft velocities up to 175 FPM, which are sufficient to disrupt the air balance of the Laminar Flow Unit.

#### E. Utilize Unidirectional Airflow - The operator must keep two important facts in mind:

- 1) The air, as supplied to the work area through the HEPA filter is contaminant-free.
- 2) Airborne contamination generated in the work area is controlled by the unidirectional flow of parallel air streams.

A solid object placed in a laminar air stream will disrupt the parallel flow and consequently, the capability of controlling lateral movement of airborne particulates. A cone of turbulence extends behind the object and laminarity of the air stream is not regained until a point is reached downstream, approximately equal to three to six times the diameter of the object. Within the parameters of this cone, particles may be carried laterally by multidirectional eddy currents.

## 5.3 Ergonomics

Ergonomics, the study or accommodation of work practices is extremely important for proper workstation usage and user health and safety. An evaluation of normal work practices should be performed with each user when working in a workstation. Evaluation criteria should be at a minimum:

- a. Proper user posture
- b. Effective workzone layout for work practice
- c. Vision or sightlines

For each of the above evaluation criterion, several aids may be supplied to accommodate the user.

- Ergonomic chair - A six-way articulating seat and back control for personalized adjustment to assure proper user posture. Be sure feet are resting on the floor, chair foot support or footrest. Also be sure back is fully supported with proper chair adjustments.
- Forearm/elbow support - The workstation is provided with a forearm support on the work access opening. Periodic mini-breaks during work practice should be taken resting forearm to avoid stress and fatigue. Elbow rests are optional that can provide support for particular work practices, such as pipetting.

- Effective workzone layout - Always prepare your work procedure to minimize reach to avoid neck and shoulder stress and fatigue. Rotating tables are optional to maximum workzone and minimize reach.
- Vision and sightline - Always prepare your work procedure to eliminate glare and bright reflections on the window. Keep your window clean and sightlines clear to your effect workzone.

## 5.4 Cleaning Procedures

### 5.4.1 General

Cleaning laboratory equipment is important in terms of both functionality and general good housekeeping. The information provided below is intended to aid the development of facility Standard Operating Procedures (SOP's) for cleaning the equipment. It is strongly recommended that all cleaning materials used be tested and verified in terms of both effectiveness and material compatibility before they are written into the cleaning SOP documentation.

- The airflow blower should be operating during the cleaning process to maintain sterility and/or containment during the cleaning process.
- Avoid directly spraying the diffuser screen while installed as this may damage the HEPA filter. Remove the diffuser screen and clean separately outside the unit.
- Apply appropriate cleaning material or surface disinfectant to surfaces. Most surface disinfectants require a specific contact time depending the materials used within the work zone. **CONSULT APPROPRIATE DISINFECTANT DOCUMENTATION FOR PROPER APPLICATION AND SAFETY PRECAUTIONS.**
  - Polycarbonate (Covestro® Makrolon®AR) has noted material compatibility concerns (see polycarbonate compatibility section). They recommend the use of Hydrogen Peroxide based materials such as the following:
    - Steriplex SD
    - Safetec surface wipes
    - Peridox RTU

It is recommended to AVOID the use of cleaning materials containing Chlorine, Quaternary Ammoniums and Phenols.

If the polycarbonate is lightly scratched, it may be able to be polished out with Mirror Glaze Plastic Polish or similar.

Further information may be available from [www.covestro.com](http://www.covestro.com)

- Stainless steel (type 304) has noted material compatibility concerns with Acids, Chlorides and Halogens. **IF THESE MATERIALS ARE USED AND ALLOWED TO BE LEFT ON THE STAINLESS STEEL SURFACE, OXIDATION AND DEGRADATION WILL OCCUR.** Only by re-wiping surfaces with either sterile water or 70% IPA will remove harmful materials from the stainless-steel surface.


Further information is available at the following:

[http://www.parrinst.com/wp-content/uploads/downloads/2011/07/Parr\\_Stainless-Steels-Corrosion-Info.pdf](http://www.parrinst.com/wp-content/uploads/downloads/2011/07/Parr_Stainless-Steels-Corrosion-Info.pdf)

**NOTE: NuAire does not offer any product warranty with respect to cleaning material compatibility. USE AT YOUR OWN RISK! The information provided above is from raw material suppliers and known general source documents for use to develop application cleaning SOP's.**

**NOTE: When cleaning the work area for the first several times, the new metal surfaces may produce some dark discoloration on the white cleaning wipes. Repeated cleaning will continuously reduce the amount of the discoloration material on the cleaning wipes over time.**

## 6.0 General Maintenance

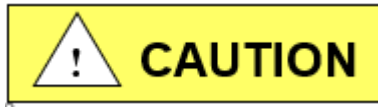

All maintenance actions on this equipment must be performed by a qualified technician who is familiar with the proper maintenance procedures required for this equipment. This includes both certification as well as repair.

### 6.1 LED/Fluorescent Lamp Replacement

The LED bulbs or older style fluorescent lamps are T8, cool white, and electronic start. The life rating is 9000 hours for a fluorescent lamp and 50,000 hours for an LED bulb based on three hour burning cycles. The method for replacing both LED and fluorescent bulbs is the same.


To replace a lamp:

1. Make sure the workstation light is off.
2. Remove (1) screw at each side of the top panel. Pull out the panel and prop it up to access the lamps.
3. Lamps are removed by rotating until the pins can be pulled down.
4. Reverse the procedure to reinstall lamps and reassemble unit.

### 6.2 HEPA Filter Replacement

The HEPA filter, under normal usage and barring an accident (puncture), does not need replacement until the efflux velocity cannot be maintained at 90 LFPM (.457 m/s)  $\pm$  10%. This may permit the HEPA filter efflux average to be as low as 81 LFPM (.411 m/s), as long as no point falls below 70 LFPM (.355 m/s). Use only replacement filters of the same rated flow and size as originally installed to ensure proper airflow and HEPA filter lifetime can be achieved.

The HEPA Filter replacement procedure is performed by the following steps: Drawing BCD-19827):

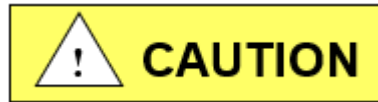

**Disconnect electrical power from the cabinet before attempting any maintenance action.**

1. Remove IV bars if installed.
2. Remove the unit diffuser screen(s) by lifting slightly and tipping out the diffuser screen(s) to remove through the front of the unit.
3. Remove filter clamp by removing the #8 x 2.00" screw, #8 washers and compression springs (Detail "A" 4 places).
4. Slightly lift and tip out filter and remove through front of the unit.

**NOTE: This step usually requires 2 people for the removal and replacement of the HEPA filter.**

5. Replace HEPA filter by tipping the filter into place resting on plenum bottom.
6. Reattach filter brackets and hardware 4 places and tighten uniformly.
7. Replace the unit diffuser screen(s) by tipping into place, lifting slightly to all screen bottom to fall into groove in work surface and to fit screen to catch the top angle to rest in place.

### 6.3 Pre-Filter Replacement

The replacement interval depends on the contaminant (large particles or lint) in the room -- a typical period is every 6 months. The prefilters are located on the top of the workstation.


### 6.4 Airflow Control System Setup and Calibration

#### 6.4.1 General


The operation of the NU-240 requires that the setup and calibration procedures be performed in order to certify or commission for usage. The setup and calibration procedures performed **ONLY BY THE CERTIFIER** ensure that setpoints are verified and that the airflow monitor sensor is calibrated to the correct values.

## 6.4.2 Configuration Parameters

Configuration parameters identify NU-240 motor type and size for proper performance characteristics. **Bold items represent default, unconfigured parameters. These are setup for the NU-240/E models at the factory, but will return to default if the control board is reset.**

- Select/Verify model and size
    - Press and hold  key, then press blower – blower – blower keys sequentially. LED indicators above both blower and UV light keys will blink fast.
- Review model type/size and change if desired as described below

- Press  $\uparrow$  or  $\downarrow$  keys to scroll through the LED segment choices associated with model type/size.

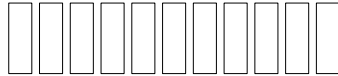


Class II, 3 foot



**Class II, 4 foot**



Class II, 5 foot

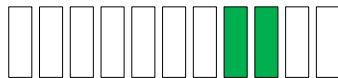


Class II, 6 foot




Laminar, 3/4 foot (NU-240/E models)




Laminar, 5/6 foot (NU-240/E models)

- Press hidden key to accept model/size and exit
  - Press  key at any time to abort and exit
- Select motor control function/type

- Press and hold  key, then press light – light – light keys sequentially. LED indicators above both blower (red) and light keys will blink fast.
- Review motor control function/type and change it desired as described below.
- Press  $\uparrow$  or  $\downarrow$  keys to scroll through the LED segment choices associated with motor control function/type.




**Auto DC ECM motor**

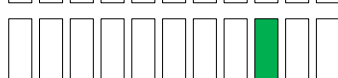

Fixed duty DC ECM motor (NU-240/E models)



Fixed duty EC Impeller motor

- Press hidden key to accept motor control function/type and exit
- Press  key at any time to abort and exit.



### 6.4.3 Airflow Calibration



Failure to calibrate airflow to the specified requirements may result in unsafe conditions of performance (i.e. product protection, noise and vibration)


The NU-240 Airflow Calibration consists of adjusting the airflow.

**THIS WORK SHOULD BE DONE ONLY BY A QUALIFIED TECHNICIAN WHO CAN MEASURE THE AIRFLOW WITH A SUITABLE VELOMETER.** NuAire provides one adjustment to set the airflow within the cabinet.

This is: PWM signal adjust via calibration parameter menu.

The PWM signal or blower speed adjustment establishes the NU-240 total volume of airflow as well as makes up for filter resistance tolerances.

DC ECM motor PWM signal DC voltage should also be monitored and recorded upon final calibration. The DC voltage may be measured using a digital voltmeter. The two test points to measure DC ECM motor voltage are located on the DC motor connector on the main control board.

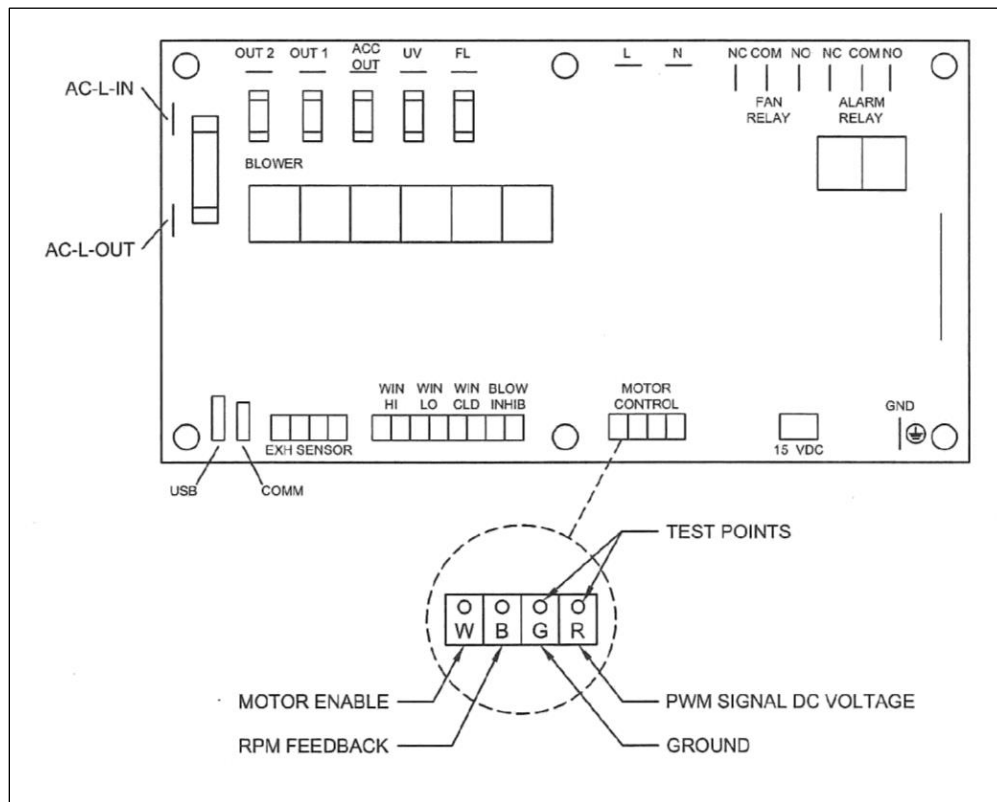

The NU-240 is considered to be certifiable if the following airflow measurements are present:


Inflow average: 90 LFPM  $\pm$  10 LFPM (.46 m/s  $\pm$  0.05 m/s) using the direct inflow measurement method or related value using the calculated inflow velocity measurement method.

**BEFORE STARTING AIRFLOW CALIBRATION PROCEDURE, LET THE NU-240 RUN FOR AT LEAST 5 MINUTES.**

### 6.4.3.1 Airflow Calibration

Step 1: Measure the inflow velocity using the recommended procedure found in 6.4.3.4. If necessary, adjust to achieve the correct average miniflow velocity within the stated range of  $90 \pm 10$  LFPM ( $.46 \pm 0.05$  m/s).

Step 2: If necessary, enter active blower speed adjustment.

- Press and hold  key, then press hidden – blower -  $\uparrow$  keys sequentially.

LED indicator above blower (green) key will blink fast.

Step 3: Press  $\uparrow$  or  $\downarrow$  keys to adjust blower speed.

- LED segments will indicate blower speed percentage and active blower speed adjustment.




- Right end red LED indicates active blower speed adjust  
The red LED will blink as soon as any adjustments are made and will continue to blink as the motor rpm settles. Once the red LED stops blinking, the motor will run steady state at the new percentage.

**Note: The red LED must be non-blinking to save or exit.**

- Green LED's indicate percentage on of scale (0-100%)
- Yellow LED's indicate minimum (left/maximum (right) blower speed has been achieved

**Note: At any time during the process:**

- Press hidden key to accept and enter the blower speed calibration point  
(If the blower speed calibration point was not successfully entered; a half second audible alarm will occur. The calibration process must then be repeated for successful entry of blower speed calibration point.)
- Press  key to abort and exit


### 6.4.3.2 PresurFlow™ Alarm Set Points

The PresurFlow™ alarm setpoints are based on the calibration setpoint. Once the calibration setpoint is entered, based on a nominal inflow velocity of 90fpm (.46mls) the associated pressure sensor value is entered as the nominal pressure value. **The high and low alarm setpoints are factory verified and set if needed at 90LFPM (.46 m/s).**

However, if specific use alarm setpoints are desired, the alarm setpoints may be adjusted by performing the following:

#### Low Alarm Setpoint


- Press and hold the Hidden and  $\downarrow$  key for 3 seconds.  
(The left red LED will blink and the green LED's indicate blower speed)
- Press  $\uparrow$  or  $\downarrow$  keys to adjust blower speed to the desired airflow velocity low alarm setpoint value.

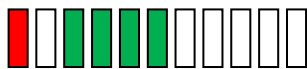

- Left end red LED indicates active low limit blower speed adjust  
The red LED will blink as soon as any adjustments are made and will continue to blink as the motor rpm settles. Once the red LED stops blinking, the motor will run steady state at the new percentage. If the low alarm setpoint value is not within an acceptable range, the left end red LED will blink at a very fast rate.

**Note: The red LED must be non-blinking to save or exit.**

**Note: At any time during the process:**

- Press hidden key to accept low alarm setpoint value
- Press and hold outlet key for three seconds to remove any previous offsets
- Press  key to abort and exit
- Upon exiting, the blower will go back to actual airflows.

### High Alarm Setpoint


- Press and hold the Hidden and ↑ key for 3 seconds.  
(The right red LED will blink and the green LED's indicate blower speed)
- Press ↑ or ↓ keys to adjust blower speed to the desired airflow velocity high alarm setpoint value.



- Red end Red LED indicates active high limit blower speed adjust.  
The red LED will blink as soon as any adjustments are made and will continue to blink as the motor rpm settles. Once the red LED stops blinking, the motor will run steady state at the new percentage.  
If the high alarm setpoint value is not within an acceptable range, the right end red LED will blink at a very fast rate.

**Note: The red LED must be non-blinking to save or exit.**



**Note: At any time during the process:**

- Press hidden key to accept low alarm setpoint value
  - Press and hold outlet key for three seconds to remove any previous offsets
  - Press  key to abort and exit
- Upon exiting, the blower will go back to actual airflows.


**Note: Specific use alarm setpoints or the offset pressure value from the nominal calibration point will be maintained with a new nominal calibration value. It is not necessary to re-enter the alarm setpoints after a nominal calibration.**

#### 6.4.3.3 PresurFlow Alarm Verification

The PresurFlow Alarm setpoints are based on the calibration setpoint. Once the calibration setpoint is entered, the Alarm setpoint offset pressure values will align from the calibration pressure value. The high or low alarm setpoint can be verified by measuring inflow volume/velocity while adjusting blower up or down within the Alarm Verification menu.

- Press and hold  key, then press ↑ ↓ ↑ sequentially releasing the  key after the 3 key sequence.
- Note: If blower was off while entering into the Alarm Verification Menu, the low alarm limit will immediately activate. Turn on blower; once airflow is above the low alarm limit, the alarm will turn off.
- LED segments will indicate blower speed percentage.



- Press ↑ or ↓ key to raise or lower blower speed. Alarm is active so yellow and red LED's will activate if pressure reaches the low or high alarm limit.
- Press  to exit (blower should turn off and not go through normal start up procedure).

**6.4.3.4 Airflow Velocity Profile**

A. Instruments: TSI 8355 Thermo anemometer

B. Procedure:

Air velocity readings are taken on a 12 inch (305mm) grid, in a plane parallel to and 6 inches (152mm) from the filter surface. No reading should be taken closer than 6 inches (152mm) from the inner edge of the filter frame.

C. Test Data:

|  |  |  |  |  |  |  |  |  |  |  |  |  |  |
|--|--|--|--|--|--|--|--|--|--|--|--|--|--|
|  |  |  |  |  |  |  |  |  |  |  |  |  |  |
|  |  |  |  |  |  |  |  |  |  |  |  |  |  |
|  |  |  |  |  |  |  |  |  |  |  |  |  |  |
|  |  |  |  |  |  |  |  |  |  |  |  |  |  |
|  |  |  |  |  |  |  |  |  |  |  |  |  |  |
|  |  |  |  |  |  |  |  |  |  |  |  |  |  |
|  |  |  |  |  |  |  |  |  |  |  |  |  |  |

|                     |                   |              |
|---------------------|-------------------|--------------|
| Number of Readings: | Average Velocity: | ft./min. m/s |
|---------------------|-------------------|--------------|

D. Acceptance Criteria:

1. Average Velocity = **80 to 100 fpm (.41 to .51 m/s)**
2. Individual Readings must be within  $\pm 20$  percent of the average velocity  
\_\_\_\_\_ to \_\_\_\_\_ fpm ( \_\_\_\_\_ to \_\_\_\_\_ m/s).

## 6.5 HEPA Filter Leak Test

In order to check filter and filter seal integrity, the HEPA filter media and seals must be directly accessible, by the measuring instrument. The challenge material (i.e. PAO) should be supplied in the top panel of the workstation. An upstream challenge may be taken by using the minihelic gauge pressure supply tube if desired.

The diffuser is secured by thumb screws on the sides.

Once removed the diffuser can be lifted up slightly using the small handles and removed from the cabinet.

| <b>Model Size</b> | <b>Filter Area (ft<sup>2</sup>)(m<sup>2</sup>)</b> |
|-------------------|--|
| 330               | 7.5 (.213)   |
| 336               | 9 (.255)   |
| 430               | 10 (.283)  |
| 436               | 12 (.340)  |
| 530               | 12.5 (.354)  |
| 536               | 15 (.425)  |
| 630               | 15 (.425)  |

### Laskin Nozzle Concentration Formula

$$\frac{\# \text{ Nozzles} \times 135 \text{ CFM} \times 100 \mu\text{g/L}}{\text{CFM}} = \text{Challenge Concentration } (\mu\text{g/L})$$

$$\frac{\# \text{ Nozzles} \times 229 \text{ CMH} \times 100 \mu\text{g/L}}{\text{CMH}} = \text{Challenge Concentration } (\mu\text{g/L})$$

## 6.6 Cleanliness Classification Test for Pharmacy Application

If this cabinet is going to be used within pharmacy, per USP797<sup>1</sup>, the cabinet must be tested to assure compliance to ISO 14644-1:2015, Cleanrooms and Associated Controlled Environments, Part 1: Classification of Air Cleanliness<sup>2</sup>. The cleanliness classification test is performed using a particle counter to measure particle counts within the cabinet workzone. Turn on cabinet and let warm up for several minutes. Turn on particle counter and flush out sample tubing line to remove latent particles. Set the particle counter to measure 0.5 micron or larger particles at the appropriate measuring rate.

### **“Operational Particle Count Test”<sup>3</sup>**

Position the particle counter isokinetic probe at a point 6 inches (152mm) upstream of the aseptic manipulation area (hand convergence point) and mounted so as not to interfere with the operator’s hand movement. The pharmacy operator will simulate IV manipulation during the particle count test using non-hazardous materials. A minimum of three (3) 1-minute particle counts shall be sampled and recorded while the user simulates aseptic compounding manipulations.

### **“At Rest Particle Count Test”**

Take 5 test points in 1-minute intervals on a grid, in a horizontal plane as measured approximately 6-inches (152mm) parallel to the diffuser. The grid location is designed as the diffuser center point and each corner measured 6-inches (152mm) from the inside perimeter.

Record the 5 particle count values for each of the test points over the 1-minute sample time. All final count particle concentrations and calculated 95% upper confidence limit shall not exceed 3520 particles per cubic meter (ppcm) or (100 particles per cubic feet (ppcf).

<sup>1</sup> USP28-NF23: United States Pharmacopeial Convention, Inc., 12601 Twinbrook Parkway, Rockville, MD 20852, USA, [www.usp.org](http://www.usp.org).

<sup>2</sup> ISO 14644-1:2015 Cleanrooms and Associated Controlled Environments-Classification of Air Cleanliness, International Organization for Standardization, Case Postale 56, CH-1211 Geneve 20, Switzerland

<sup>3</sup> CAG-002-2006: CETA Compounding Isolator Testing Guide, Controlled Environment Testing Association, 1500 Sunday Drive, Suite 102, Raleigh, NC 27607, USA, [www.cetainternational.org](http://www.cetainternational.org)

## 6.7 Main Control Board Description and Replacement


To access the main control board for fuse or board replacement, remove (1) rear screw at each top side of the top panel and pull out panel to access. Now the main control board is exposed for service.

### 6.7.1 Main Control Board Replacement

The main control board consists of one Printed Circuit Board (PCB) assembly.

The PCB contains the power supply, configuration switch, sensor inputs/outputs and control inputs/outputs components and display.

### 6.7.2 Main Control Board Fuse Replacement

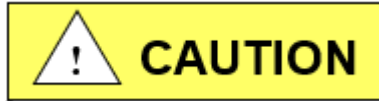


Disconnect electrical power before fuse replacement.

All AC circuits are fuse protected and when replacement is necessary, USE ONLY FUSES OF SAME TYPE AND RATING FOR PROTECTION AGAINST RISK OF FIRE.

| DESCRIPTION: | BLOWER FUSE      | OUTLET FUSE<br>(Option Only) | ACCESSORY OUTPUT FUSE<br>(Option Only) | LIGHT FUSES |
|--------------|------------------|------------------------------|--|-------------|
| FUSE TYPE:   | TIME-LAG         | TIME-LAG                     | TIME-LAG                               | TIME-LAG    |
| FUSE SIZE:   | 1/4 X 1-1/4 INCH | 5 X 20MM                     | 5 X 20MM                               | 5 X 20MM    |
| NU-240-3XX   | 10 AMPS          | 3 AMPS                       | 2 AMPS                                 | 1 AMP (2)   |
| NU-240-4XX   | 10 AMPS          | 3 AMPS                       | 2 AMPS                                 | 1 AMP (2)   |
| NU-240-5XX   | 10 AMPS          | 3 AMPS                       | 2 AMPS                                 | 1 AMP (2)   |
| NU-240-6XX   | 10 AMPS          | 3 AMPS                       | 2 AMPS                                 | 1 AMP (2)   |

### 6.7.3 Main Control Board Replacement

**Note:** All setup and calibration data will be lost, the memory reinitialized to the default values and all control functions reset to an initial cabinet power condition. If possible, before the main control board replacement, it would be preferred to know the operational parameters of the cabinet, (i.e. blower speed/PWM signal DC voltage setpoints if modified and airflow data from previous certification).

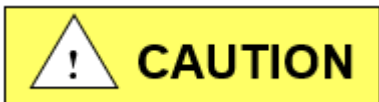

Disconnect electrical power before attempting any maintenance action.

The main control board is fastened to the control center with (6) 6-32 studs/nuts. All electrical connections are made with removable terminals and/or Faston connectors except for the motor/blower connector which uses a screw terminal. Remove all electrical connections and fasteners and then remove the main control board from the control center.


Install new main control board by reattaching all electrical connections and fasteners. Once installed, rotate control center to normal position and fasten in place. Reconnect power to cabinet.

### 6.7.4 Aeromax™ Control System Reset

The Aeromax™ control board has a reset function available for Qualified Service Personnel.

Master Reset-Resets all calibration, cabinet size, motor type/function and option settings back to their default condition. It is recommended to perform a master reset upon installation of replacement control board or if there seems to be intermittent functional abnormalities.

Perform the following sequence to accomplish a master reset:


1. Turn off all functions i.e. blower, lights, etc.
2. Press and hold alarm silence key for 10 seconds until alarm LED flashes, then release.
3. Press hidden key within 5 seconds of the alarm LED flashing. There will be a (3) beep confirmation that the reset has occurred.
4. Turn off cabinet (either unplug cabinet or use power switch on left side within control panel) then turn back on to reinitialize the Aeromax™ control board. Just as in the initial power up, the RED alarm LED will blink to indicate power up status.
5. Press  key to clear the power on reset message.
6. Enter into service parameters and input configuration for calibration and option settings as required.

## 7.0 Option Parameters

The option parameter menu allows **A QUALIFIED TECHNICIAN** to configure several different optional parameters per the menu as described below.


### 7.1 Sync Function with Active Blower


To access the option parameter menu, perform the following:

- Press and hold  key, then press hidden - Blower – LED Light keys sequentially. Red LED indicator above the blower key will blink fast

The PresurFlow™ blinking green LED segments will indicate seven optional parameters as shown and described below. The UV Light key (move left) and outlet key (move right) allows selection of the option parameter desired.

Once the desired option parameter is indicated, press  $\uparrow$  or  $\downarrow$  key to turn on or off. A slow blinking green LED indicator means off and a fast blinking green LED indicator means on. Multiple option changes can be selected.

- Pressing the hidden key will accept all changes and exit
- Pressing the  key will abort the process and exit

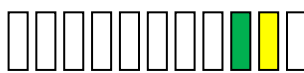

Sync Fan Relay with Active Blower -

Normally the fan relay will activate when the blower switch is pressed. Blower can either be actively running or pending. If the fan relay sync is active the blower must be actively running for the relay to change state.



Sync Accessory Outlet with Active Blower - (only usable if accessory outlet is installed)

Normally the accessory outlet is on all the time. If the accessory outlet sync is active, the blower must be actively running for the accessory outlet to turn on.



Sync Outlet Power with Active Blower - (only usable if outlet is hardwired to control board)

Normally the outlet power is turned on via the outlet key. If the outlet power sync is active, the outlet power will turn on and off with the blower or may be turned on and off independently if the blower is active.



Sync LED Light with Active Blower -

Normally the LED light is turned on via the LED light key. If the LED light sync is active, the LED light will turn on and off with the blower or may be turned on and off independently if the blower is active.




Sync 15 Volt DC output with active blower normally the 15 Volt DC output located on

the control board is on when power is applied to the system. If the 15 Volt DC output sync is active, the blower must be actively running for the 15 Volt DC output to turn on.

### 7.2 Blower/Airflow Options


To access the option parameter menu, perform the following:

- Press and hold  key, then press hidden -  $\uparrow$  and  $\downarrow$  keys sequentially. Red LED indicator above the blower key will blink fast

The PresurFlow™ blinking green LED segments will indicate seven optional parameters as shown and described below.

The UV Light key (moves left) and outlet key (moves right) allows selection of the option parameter desired.

Once the desired option parameter is indicated, press  $\uparrow$  or  $\downarrow$  key to turn on or off. A slow blinking green LED indicator means off and a fast blinking green LED indicator means on. Multiple option changes can be selected.


- Pressing the hidden key will accept all changes and exit
- Pressing the  key will abort the process and exit



**Require Password –**

Normally it is not required to use a password (i.e. 3 key press sequence of the blower and hidden key). If the option is turned on, it would be required to use the correct password to turn on the blower. The default password once turned on is blower-hidden-blower keys in sequence. The password can be changed in the blower password option menu.

Note: If turning off the password option, you must also select the no password menu item (see section 5.1.8.3)

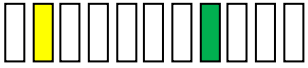

**Allow UV light anytime – (not available on NU-240 products)**


Normally the UV light is interlocked with the window being in the closed position. For service purposes only, if the function is active, the UV light may be turned on at any window height.

**Note: In addition to the Allow UV light anytime system function, there is a double redundant UV light window interlock relay. To override the UV light window interlock relay, the relay itself must also be shorted. (See electrical schematic for reference).**



**Manual Blower Restart –**

Normally when the blower is actively running and a power interruption occurs. The blower will automatically come back on when power is restored. If this function is turned on, the blower will not automatically come back after a power interruption, but would require the user to press the blower key to restart the blower.

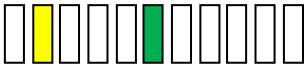

**Low Flow –**

Normally the function is turned off. If selected and turned on, once the blower is actively running. Upon closure, the blower will continue to run at a calibrated lower speed level to maintain a negative airflow. The PresurFlow™ will indicate(s) blinking green LEDs along with green LED above blower key.



**Disable PresurFlow™ -**

If this function is active, the digital pressure sensor and alarm function are turned off.




**Temporary Low Flow – (Function not available on NU-240 models.)**

When this option parameter is turned on and the other requirements below are met, the blower key (when held for 5 seconds) will toggle blower between normal and Low Flow blower speed. The Low Flow blower will time out, based on the Auto Timer duration for Low Flow without a blower key press and the blower speed will revert back to normal.

In addition to turning this option on the following requirements must also be met

- Temporary Low Flow Option must be selected
- Password Option must be selected
- Window is at normal height and blower is running
- Low Flow Option must be selected
- Low Flow blower auto timer must be set for a minimum of 15 minutes.

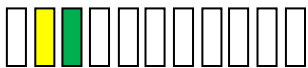

**Disable audible alarms / audible key feedback**

Normally audible alarms from the PresurFlow™ monitor and audible key feedback are present to provide an audible sound to the user. If this function is active, all audible sound will be silenced when the control system is in normal run mode. Audible sound will still occur in any service/calibration menu function.



## 8.0 Error Messages and Troubleshooting

Audible alarms and error messages occur for a variety of reasons. Whenever an alarm condition has been present for a period of at least 10 seconds, the audible alarm/error message will be presented and stay on until the error is cleared. The audible alarm will be on for 30 seconds upon initial alarm condition, then once every ten seconds. When presented with an error message, please perform the following:

Step 1: NOTE ALL ERROR MESSAGES.

Error message will appear on the control panel with red LED's.

Step 2: VERIFY ERROR MESSAGES.

Error messages can be verified by clearing the error function by either turning the blower or the cabinet on and off.

Step 3: MONITOR RE-OCCURRENCE OF ERROR MESSAGES.

If re-occurrence of the error message is immediate or daily, use the following guide to correct the situation.

| Issue   | Error Description   | Correction   |
|---|---|--|
| Cabinet LED lights won't turn on.                         | Blue LED above light key indicates the lamp should be on.                       | Check light fuse on main control board.<br>Check LED lamps.<br>Check voltage coming out of main control board to light ballasts.<br>Check light starters, if present.<br>Check ballast.              |
| Cabinet blower won't turn on.                             | Green LED above blower key indicates the blower should be on.<br>Airflow Alarm. | Check blower fuse on main control board.<br>Check AC voltage coming out of main control board.<br>Check wiring to blower.<br>Check blower motor.<br>Check DC motor PWM signal on main control board. |
| Red alarm LED blinks                                      | Indicates a power interruption has occurred.                                    | Press any key to clear.  |
| Blower or light fuse continues to blow after replacement. | N/A   | Check for short on output of fuse.<br>Isolate output of fuse by disconnecting control center connectors, light circuit, AC or DC blower circuit, etc. to isolate the short.                          |
| PresurFlow™ left red LED indicator on and red LED alarm   | PresurFlow™ reading low flow (pressure)   | Check airflow values.<br>Check blower function.<br>Recalibrate PresurFlow™ system.   |
| PresurFlow™ right red LED indicator on and red LED alarm  | PresurFlow™ reading high flow (pressure)  | Check airflow values.<br>Recalibrate PresurFlow™ system.   |
| All PresurFlow™ LED's blink                               | Message acknowledges new firmware was loaded into microprocessor                | N/A  |
| Blower red LED blinks and red LED alarm                   | Indicates that the motor rpm signal has been interrupted                        | Check connectors and wires from main control board to the motor<br>Replace motor if required   |

# Environmental resistance of Makrolon® polycarbonate sheet

**Makrolon® POLYCARBONATE SHEET IS RESISTANT AT 70°F AND 0% STRAIN TO:**

**CHEMICALS:**

|                        |                         |                       |                        |                      |
|------------------------|-------------------------|-----------------------|------------------------|----------------------|
| Amyl Alcohol           | Chromic Acid (20%)      | Lactic Acid (20%)     | Potassium Bromate      | Sodium Chloride      |
| Aluminum Chloride      | Citric Acid (40%)       | Magnesium Chloride    | Potassium Bromide      | Sodium Hypochlorite  |
| Aluminum Sulfate       | Copper Chloride         | Magnesium Sulfate     | Potassium Nitrate      | Sodium Sulfate       |
| Ammonium Chloride      | Copper Sulfate          | Manganese Sulfate     | Potassium Perchlorate  | Stannous Chloride    |
| Ammonium Nitrate       | Formic Acid (10%)       | Mercuric Chloride     | Potassium Permanganate | Sulfur               |
| Ammonium Sulfate       | Formalin (30%)          | Nickel Sulfate        | Potassium Persulfate   | Sulfuric Acid (10%)* |
| Antimony Trichloride   | Glycerine               | Nitric Acid (10%)     | Potassium Sulfate      | Sulfuric Acid (50%)  |
| Arsenic Acid           | Heptane                 | Nitric Acid (20%)     | Silicone Oil           | Tartaric Acid (30%)  |
| Butyl Alcohol          | Hydrochloric Acid (10%) | Oleic Acid            | Silver Nitrate         | Zinc Chloride        |
| Calcium Nitrate        | Hydrogen Peroxide (30%) | Oxalic acid           | Sodium Bicarbonate     | Zinc Sulfate         |
| Chlorinated Lime Paste | Hydrofluoric Acid (10%) | Pentane               | Sodium Bisulfate       |                      |
| Chrome Alum            | Isopropanol             | Phosphoric Acid (10%) | Sodium Carbonate       |                      |

\*Sulfuric Acid at 1% attacks polycarbonate sheet

**COMMON HOUSEHOLD MATERIALS Makrolon® POLYCARBONATE IS RESISTANT TO:**

**CHEMICALS:**

|                  |                      |                         |
|------------------|----------------------|-------------------------|
| Borax            | Joy Liquid Detergent | Rum                     |
| Cocoa            | Insulating Tape      | Salad Oil               |
| Cement           | Linseed Oil          | Salt Solution (10%)     |
| Chocolate        | Liquor               | Soap (Soft/Hard)        |
| Cod Liver Oil    | Milk                 | Table Vinegar           |
| Cognac           | Mineral Water        | Tincture of Iodine (5%) |
| Coffee           | Mustard              | Tomato Juice            |
| Detergents       | Olive Oil            | Vodka                   |
| Fish Oil         | Onions               | Washing Soap            |
| Fruit Syrup      | Orange Juice         | Water                   |
| Grapefruit Juice | Paraffin Oil         | Wine                    |
| Gypsum           | Rapeseed Oil         |                         |

**PETROLEUM PRODUCTS Makrolon® POLYCARBONATE SHEET IS RESISTANT TO:**

|                |                 |
|----------------|-----------------|
| Compressor Oil | Spindle Oil     |
| Diesel Oil     | Transformer Oil |
| Kerosene       | Vacuum Pump Oil |
| Refined Oil    |                 |

Note: Elevated temperature and/or strain significantly alters resistance to industrial petroleum products.

**LIMITED RESISTANCE AT 70°F AND 0% STRAIN TO:**

|                  |                           |
|------------------|---------------------------|
| Antifreeze       | Hydrochloric Acid (conc.) |
| Calcium Chloride | Milk or Lime (CaOH)       |
| Cyclohexanol     | Nitric Acid (conc.)       |
| Ethylene Glycol  | Sulfuric Acid (conc.)     |

**Makrolon® POLYCARBONATE SHEET IS NOT RESISTANT TO:**

**CHEMICALS:**

|                     |                          |                       |                                |                        |
|---------------------|--------------------------|-----------------------|--------------------------------|------------------------|
| Acetaldehyde        | Benzyl Alcohol           | Chlorobenzene         | Formic Acid (conc.)            | Phosphorus Trichloride |
| Acetic Acid (conc.) | Brake Fluid              | Chloroethene          | Freon (refrigerant/propellant) | Propionic Acid         |
| Acetone             | Bromobenzene             | Cutting Oils          | Gasoline                       | Sodium Sulfide         |
| Acrylonitrile       | Butylic Acid             | Cyclo Hexanone        | Lacquer Thinner                | Sodium Hydroxide       |
| Ammonia             | Carbon Tetrachloride     | Cyclohexene           | Methyl Alcohol                 | Sodium Nitrate         |
| Ammonium Fluoride   | Carbon Disulfide         | Dimethyl Formamide    | Nitrobenzene                   | Tetrahydronaphthalene  |
| Ammonium Hydroxide  | Carbolic Acid            | Ethane Tetrachloride  | Nitrocellulose Lacquer         | Thiophene              |
| Ammonium Sulfide    | Caustic Potash Sol. (5%) | Ethylamine            | Ozone                          | Toluene                |
| Benzene             | Caustic Soda Sol. (5%)   | Ethyl Ether           | Phenol                         | Turpentine             |
| Benzoic Acid        | Chloride                 | Ethylene Chlorohydrin | Phosphorus Hydroxy             | Xylene                 |

**Makrolon® POLYCARBONATE SHEET IS DISSOLVED BY:**

Chloroform, Cresol, Dioxane, Ethylene Dichloride, Methylene Chloride, Pyridine

**EFFECTS OF MOISTURE ON Makrolon® POLYCARBONATE SHEET:**

Makrolon Polycarbonate Sheet has good resistance to water up to approximately 150°F. Above this temperature, the effect of moisture is time-temperature related. Exposing Makrolon Polycar-

bonate Sheet to repeated steam cleaning or dishwashing can create hydraulic crazing. The result can be a clouding of the surface and ultimately a loss of physical strength properties.

## 10.0 Remote Contacts


The NU-240 has several contact closures for remote sensing of various functions.

### 10.1 Fan Relay

The fan relay contacts are normally open and closed contact closure outputs that are activated whenever the blower key is pressed and the blower key LED indicator is on or blinking. Contact ratings are 250 VAC maximum at 2 Amps.

### 10.2 Alarm Relay

The alarm relay contacts are normally open and closed contact closure outputs which are activated whenever an airflow alarm condition occurs. Contact ratings are 250 VAC maximum at 2 Amps.

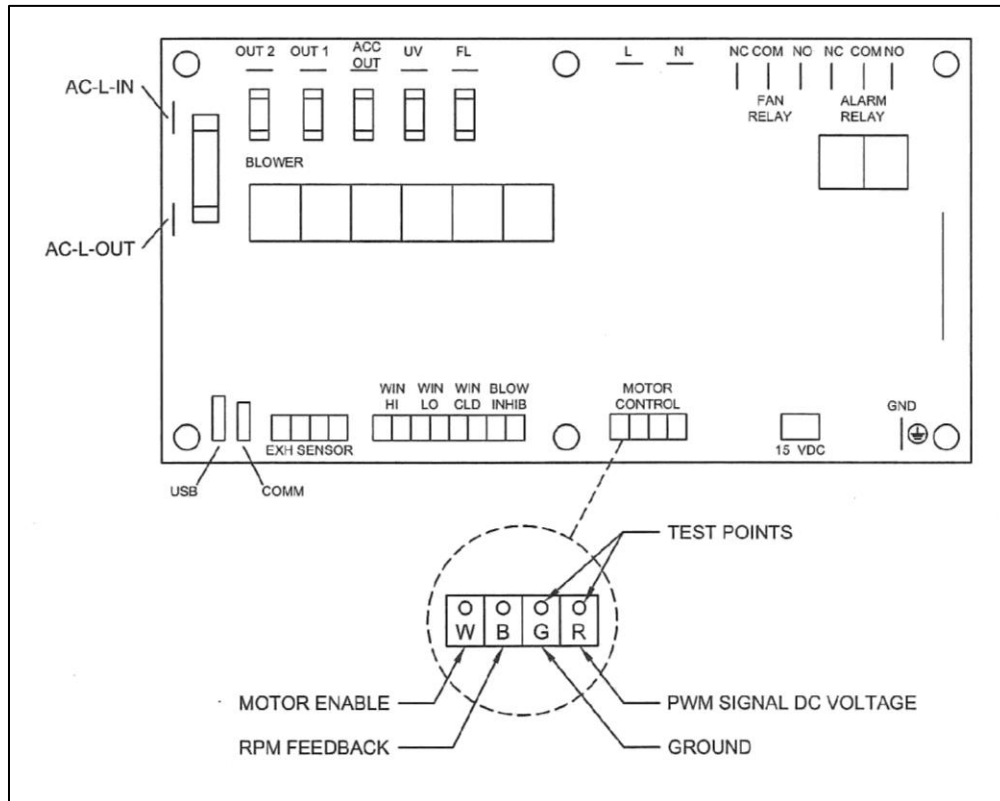

### 10.3 15VDC Output

The 15VDC (100mA) output is generated if the blower is actively running.

## 11.0 Electrical/Environmental Requirements

### 11.1 Electrical

(Supply voltage fluctuations not to exceed +/- 10%)

| <u>Electrical</u> | <u>Voltage</u> | <u>Phase</u> | <u>Frequency</u> | <u>Current (Amps)</u> | <u>Independent Outlet Current for Each Outlet (Optional)</u> |
|-------------------|----------------|--------------|------------------|-----------------------|--|
| NU-240-330        | 115 Volts      | 1            | 60 Hz            | 7                     | 10   |
| NU-240-336        | 115 Volts      | 1            | 60 Hz            | 7                     | 10   |
| NU-240-430        | 115 Volts      | 1            | 60 Hz            | 8                     | 10   |
| NU-240-436        | 115 Volts      | 1            | 60 Hz            | 8                     | 10   |
| NU-240-530        | 115 Volts      | 1            | 60 Hz            | 8                     | 10   |
| NU-240-536        | 115 Volts      | 1            | 60 Hz            | 8                     | 10   |
| NU-240-630        | 115 Volts      | 1            | 60 Hz            | 8                     | 10   |

### 11.2 Operational Performance (for indoor use only)

Environment Temperature Range: 60°F - 90°F (15.6°C - 32.2°C)

Environment Humidity: Maximum relative humidity 80% for temperatures up to 31°C decreasing linearly to 50% relative humidity at 40°C

Environment Altitude: 6562 Feet (2000 Meters) maximum

### 11.3 Light Exposure

Standard LED Lighting @ 150 ft. candles (1614 LUX) maximum intensity.

### 11.4 Installation Category: 2.0

Installation category (overvoltage category) defines the level of transient overvoltage which the instrument is designed to withstand safely. It depends on the nature of the electricity supply and its overvoltage protection means. For example, in CAT II, which is the category used for instruments in installations supplied from a supply comparable to public mains such as hospital and research laboratories and most industrial laboratories, the expected transient overvoltage is 2500V for a 230V supply and 1500V for a 120V supply. Main supply fluctuations are not to exceed  $\pm 10\%$  of nominal voltage.

### 11.5 Pollution Degree: 2.0

Pollution degree describes the amount of conductive pollution present in the operating environment. Pollution degree 2 assumes that normally only non-conductive pollution such as dust occurs with the exception of occasional conductivity caused by condensation.

### 11.6 Chemical Exposure

Chemical exposure should be limited to antibacterial materials used for cleaning and disinfecting. Chlorinated and Halogen materials are not recommended for use on stainless steel surfaces.

### 11.7 EMC Performance (classified for light industrial)

Emissions: EN61326-1 (Class A)

Immunity: EN61326-1 (Industrial electromagnetic environment)



The EMC performance requirements are generated within the product enclosure. The enclosure will be all metal grounded to earth. In addition, the membrane front panel will also include a ground plane for maximum protection and an electrostatic shield.


## 12.0 Disposal and Recycling

Workstations that are no longer in use and are ready for disposal contain reusable materials. ALL components may be disposed and/or recycled after they are known to be properly disinfected.

**NOTE: Follow all local, state, and federal guidelines for disposal of HEPA filter solid waste.**



Recycle

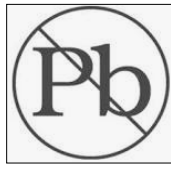

Lead Free

| <b>Component</b>        | <b>Material</b>                 |
|-------------------------|---------------------------------|
| Main Workstation        | Painted Steel                   |
| Worksurface             | PVC / Stainless Steel           |
| Side Panels             | Polycarbonate / Stainless Steel |
| HEPA Filter Frame       | Aluminum                        |
| Control Center          | Painted Steel                   |
| Diffuser                | Aluminum                        |
| HEPA Filter             | Aluminum                        |
| Hepex Bag               | PVC                             |
| Blower Wheel & Housing  | Steel                           |
| Motor                   | Various Steel / Copper          |
| Printed Wiring Assembly | Lead Free Electronic            |
| Wire                    | PVC Coated Copper               |
| Ballasts                | Various Steel, Electronic       |
| Connectors              | Nylon                           |
| Hardware                | Stainless Steel and Steel       |


**Note: Material type can be verified with use of a magnet with stainless and aluminum being non-magnetic.**

| REV | CO     | DESCRIPTION          | DATE       | DFTR | CHKD |
|-----|--------|----------------------|------------|------|------|
| C   | 002148 | NEW FRONT PANEL ASSY | 05/15/2023 | JB   | LS   |

HEPA FILTER REPLACEMENT PROCEDURE

1. DISCONNECT UNIT POWER.
2. REMOVE IV BARS IF INSTALLED.
3. REMOVE THE UNIT DIFFUSER SCREEN(S) BY LIFTING SLIGHTLY AND TIPPING OUT THE DIFFUSER CLAMP BY REMOVING THROUGH THE FRONT OF THE UNIT.
4. REMOVE FILTER CLAMP BY REMOVING THE #6 X 2.00" SCREW, #8 WASHERS AND COMPRESSION SPRINGS (DETAIL "A" 4 PLACES).
5. SLIGHTLY LIFT AND TIP OUT HEPA FILTER AND REMOVE THROUGH FRONT OF UNIT.
6. REPLACE HEPA FILTER BY TIPPING THE FILTER INTO PLACE RESTING ON PLENUM BOTTOM.
7. REATTACH FILTER BRACKETS AND HARDWARE 4 PLACES AND TIGHTEN UNIFORMLY.
8. REPLACE THE UNIT DIFFUSER SCREEN(S) BY TIPPING INTO PLACE, LIFTING SLIGHTLY TO ALLOW SCREEN BOTTOM TO FALL INTO GROVE IN WORK SURFACE AND TO FIT SCREEN TO CATCH THE TOP ANGLE TO REST IN PLACE.
9. RECONNECT UNIT POWER.



| ITEM | QTY | PART NUMBER  | DESCRIPTION   |
|------|-----|--------------|---|
| 1    | 2   | A-980356-17  | Prefiller 16X20 (Black 10PPH) Foam                  |
| 2    | 1   | A-980979     | 36X46X3 Aluminum Frame HEPA Filter Bead Pack, 99.99 |
| 3    | 1   | B-240-999    | DIFFUSER PANEL                                      |
| 4    | 4   | B-240-1005   | Filter Clamp Bracket                                |
| 5    | 4   | X-980267-03  | Die Spring 1/2 X 1 Medium Heavy Duty                |
| 6    | 8   | X-980794-01  | Washer #8 SST                                       |
| 7    | 4   | X-980848-03  | Screw HD Socket #6-32 X 2 SST                       |
| 56   | 4   | X-980893-02  | Weldpin 1/4 X 1/2                                   |
| 57   | 1   | AF-000220-02 | Front Panel Assembly 4'                             |
| 58   | 1   | AF-000224-02 | Top Frame Assembly 4'                               |

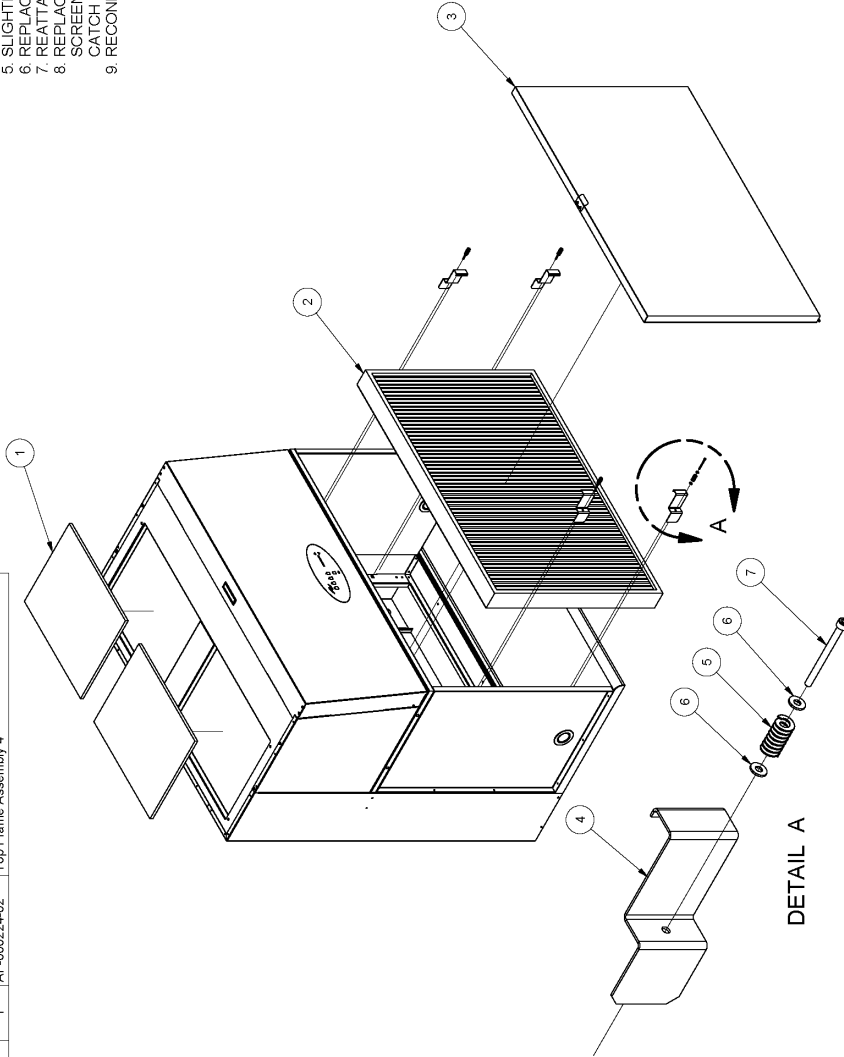

|   |           |
|---|-----------|
|   |           |
| TITLE NU-240 HORIZONTAL AIRFLOW FILTER REPLACEMENT INSTRUCTIONS |           |
| DTM   | DATE      |
| LS  | 3/20/2018 |
| DRAWING NUMBER  | BCD-19827 |
| CRD SHEET   | 1 OF 1    |
| BP  |           |


| REV | ECO    | DESCRIPTION                  | DATE       | DFTM | CHKD |
|-----|--------|------------------------------|------------|------|------|
| B   | 000529 | UPDATED TITLE AND HOLE NOTES | 12/15/2020 | LS   | BP   |

| PARTS LIST |     | DESCRIPTION |
|------------|-----|-------------|
| ITEM       | QTY | PART NUMBER |
| 1          | 2   | X-980610-01 |
| 2          | 2   | X-980794-02 |
| 56         | 4   | X-980803-02 |

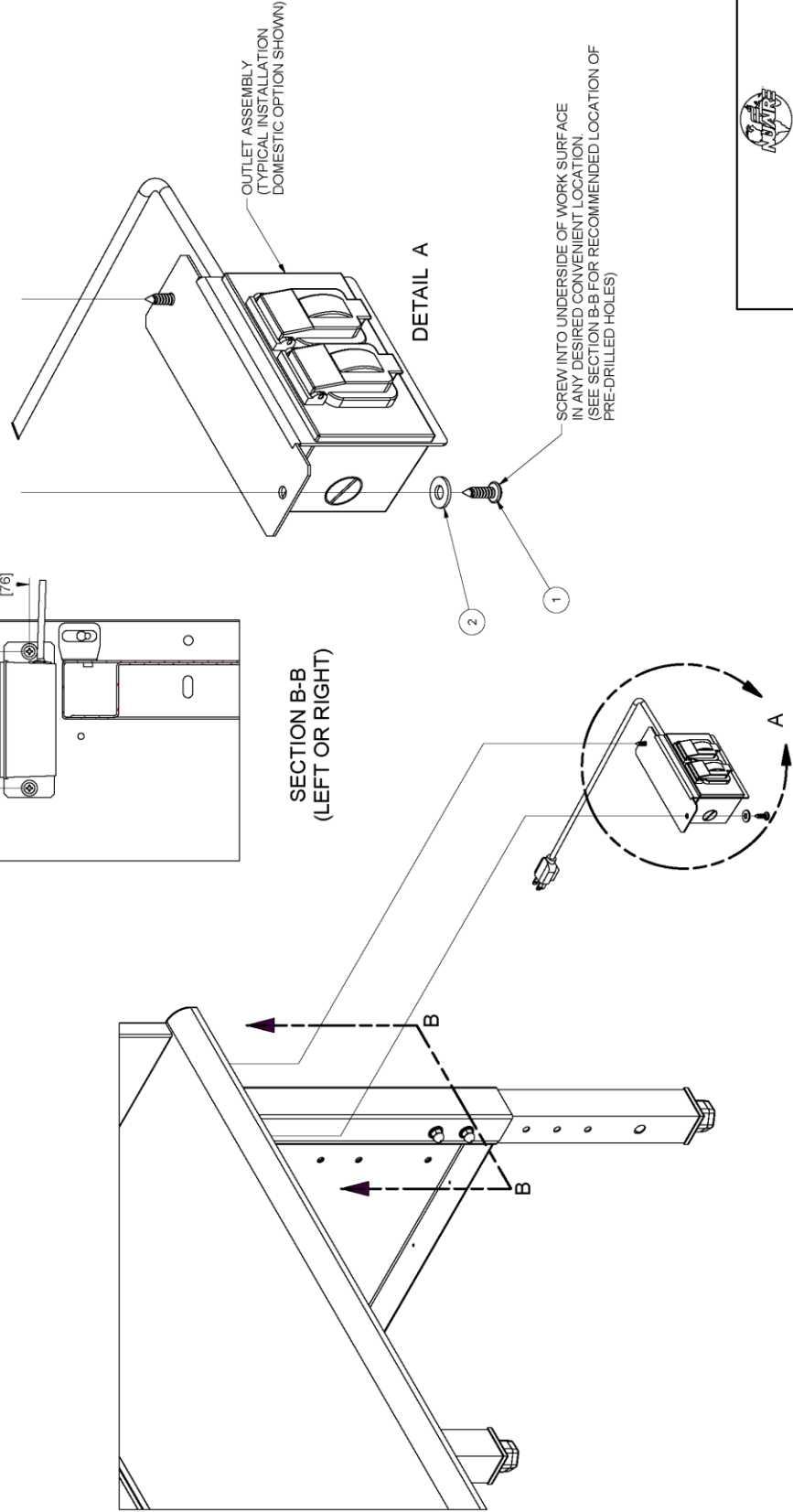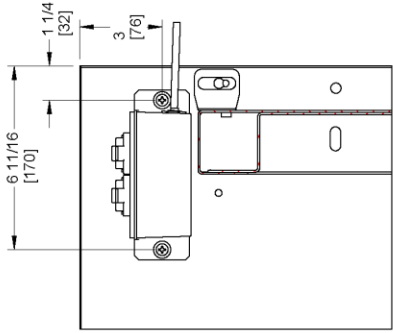

|   |                |
|---|----------------|
|   |                |
| TITLE NU-240/340 Outlet Installation Instructions |                |
| DFTM TH   | DATE 9/12/2018 |
| CRKD  | SHEET 1 OF 1   |
| DRAWING NUMBER                                    | BCD-20293      |
|   | B              |

| REV | CO     | DESCRIPTION         | DATE      | DFTR | CHKD |
|-----|--------|---------------------|-----------|------|------|
| E   | 002352 | UPDATED MOTOR FUSES | 8/11/2023 | LS   | JS   |

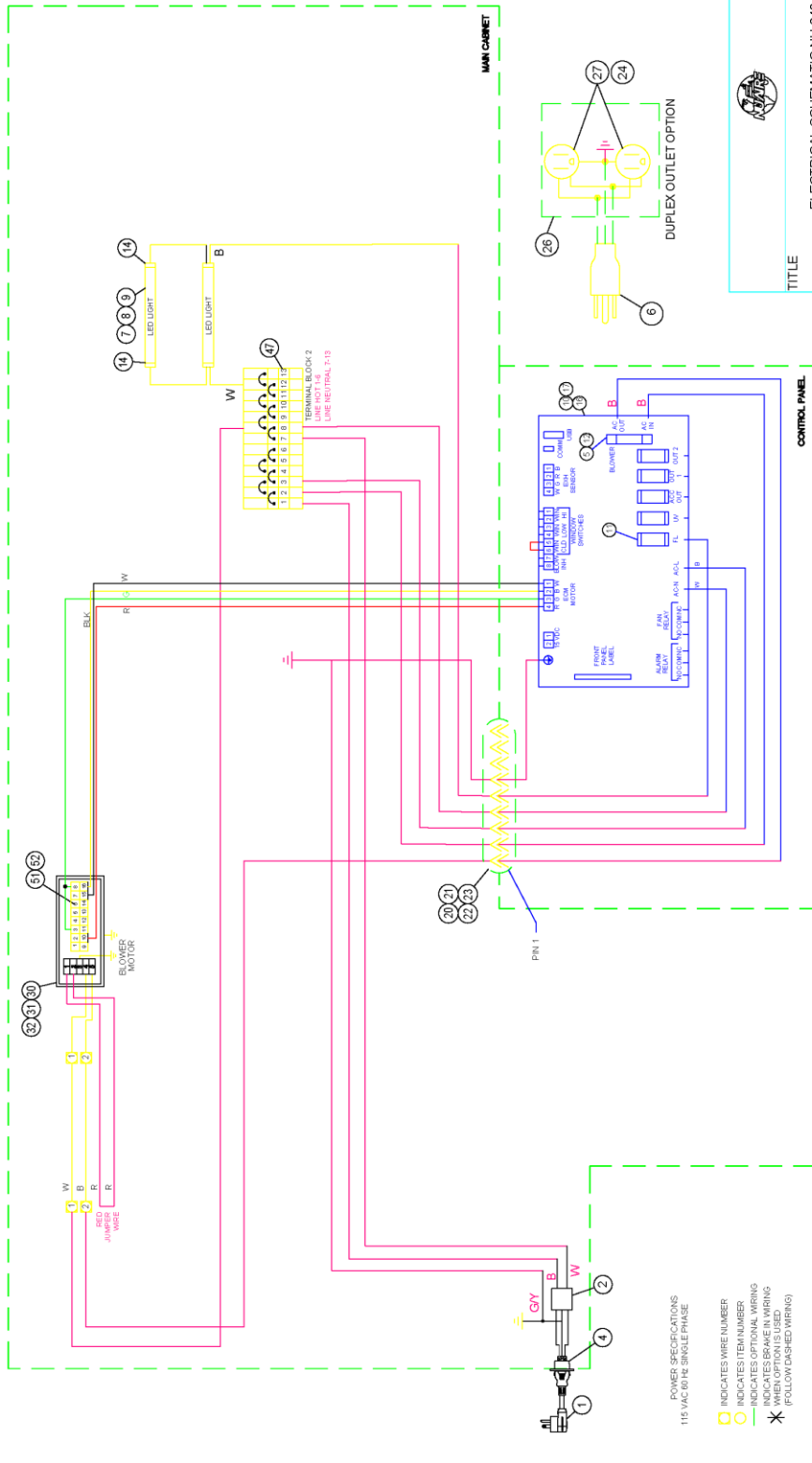


TITLE  
ELECTRICAL SCHEMATIC NU-240 115V 60 HZ  
(SERIES B)

| DFTR | TH | DATE      | CHKD | BK | SHEET | 1 OF 2 |
|------|----|-----------|------|----|-------|--------|
|      |    | 3/23/2021 |      |    |       |        |

DRAWING NUMBER  
CD-000098



| ITEM | QTY | DESCRIPTION   | NVAIRE P/N   | REV | CO     | DESCRIPTION | DATE | DFTR | CHKD |
|------|-----|---|--|-----|--------|-------------|------|------|------|
| 1    | 1   | POWER CORD, NE/IA 5-15P, IEC C18 AUTO-LOCK  | 996-00018  | E   | 002352 | SEE SHEET 1 |      |      |      |
| 2    | 1   | EMI FILTER  | X-599319-03  |     |        |             |      |      |      |
| 3    |     |   |  |     |        |             |      |      |      |
| 4    | 1   | POWER CORD CONNECTOR, IEC C28   | 996-00017  |     |        |             |      |      |      |
| 5    | 1   | FUSE (10A) (14X1-1/4) (MOTOR ALL SIZES)   | X-999921-01  |     |        |             |      |      |      |
| 6    | A/R | 15A POWER CORD 12 FT. INULDED PLUG  | X-999543   |     |        |             |      |      |      |
| 7    | A/R | LED BYPASS BULB 48"   | X-999914-01  |     |        |             |      |      |      |
| 8    | A/R | LED BYPASS BULB 48"   | 996-00014-02   |     |        |             |      |      |      |
| 9    | A/R | LED BYPASS BULB 48"   | 996-00014-03   |     |        |             |      |      |      |
| 10   | 1   | ASSY, GAF III (BASIC) (3 FT.)   | F-340754   |     |        |             |      |      |      |
| 11   | 1   | FUSE (1 AMP) (5 X 20mm) (LIGHTS)  | X-999035-13  |     |        |             |      |      |      |
| 12   |     |   |  |     |        |             |      |      |      |
| 13   |     |   |  |     |        |             |      |      |      |
| 14   | 4   | LAMPHOLDER FIXED BASE B/PIN   | X-999163-02  |     |        |             |      |      |      |
| 15   |     |   |  |     |        |             |      |      |      |
| 16   | 1   | ASSY, GAF III (BASIC) (4, 5 & 6 FT.)  | F-340755   |     |        |             |      |      |      |
| 17   | 1   | GAF III ASSY (BASIC; ALL SIZES)   | BEC-000542   |     |        |             |      |      |      |
| 18   |     |   |  |     |        |             |      |      |      |
| 19   |     |   |  |     |        |             |      |      |      |
| 20   | 1   | CONNECTOR - RECEPTACLE 9 PIN  | X-999234-02  |     |        |             |      |      |      |
| 21   | 1   | CONNECTOR - PLUG 9 PIN  | X-999235-02  |     |        |             |      |      |      |
| 22   | A/R | SOCKET CONTACT  | X-999244-02  |     |        |             |      |      |      |
| 23   | A/R | PIN CONTACT   | X-999244-01  |     |        |             |      |      |      |
| 24   | A/R | 15 AMP GFI OUTLET   | X-999330   |     |        |             |      |      |      |
| 25   |     |   |  |     |        |             |      |      |      |
| 26   | A/R | ELECTRICAL JUNCTION BOX   | X-999373   |     |        |             |      |      |      |
| 27   | A/R | OUTLET  | X-999511   |     |        |             |      |      |      |
| 28   |     |   |  |     |        |             |      |      |      |
| 29   | A/R | BLOWER ASSEMBLY (SERIES 1) 240-530  | AF-000196-01   |     |        |             |      |      |      |
| 30   | A/R | BLOWER ASSEMBLY (SERIES 1) 240-5364*  | AF-000196-02   |     |        |             |      |      |      |
| 31   | A/R | BLOWER ASSEMBLY (SERIES 1) 240-5364*  | AF-000196-03   |     |        |             |      |      |      |
| 32   |     |   |  |     |        |             |      |      |      |
| 33   |     |   |  |     |        |             |      |      |      |
| 34   |     |   |  |     |        |             |      |      |      |
| 35   |     |   |  |     |        |             |      |      |      |
| 36   |     |   |  |     |        |             |      |      |      |
| 37   |     |   |  |     |        |             |      |      |      |
| 38   |     |   |  |     |        |             |      |      |      |
| 39   |     |   |  |     |        |             |      |      |      |
| 40   |     |   |  |     |        |             |      |      |      |
| 41   |     |   |  |     |        |             |      |      |      |
| 42   |     |   |  |     |        |             |      |      |      |
| 43   |     |   |  |     |        |             |      |      |      |
| 44   |     |   |  |     |        |             |      |      |      |
| 45   |     |   |  |     |        |             |      |      |      |
| 46   |     |   |  |     |        |             |      |      |      |
| 47   | 1   | TERMINAL BLOCK ASSY CONTAINING:<br>• END NO PLACE, TERMINAL BLOCK<br>• END, 4 PLACE, TERMINAL BLOCK<br>• TERMINAL BLOCK, 4 PLACE NON-END<br>• TERMINAL BLOCK, 4 PLACE NON-END<br>• JUMPER, TERMINAL BLOCK<br>• JUMPER, TERMINAL BLOCK | A-435-4099-02<br>X-999921-01<br>X-999921-03<br>X-999919-01<br>X-999920-01<br>X-999920-02<br>X-999920-01<br>X-999920-01 |     |        |             |      |      |      |
| 48   | 3   | JUMPER, TERMINAL BLOCK  | X-999920-01  |     |        |             |      |      |      |
| 49   |     |   |  |     |        |             |      |      |      |
| 50   |     |   |  |     |        |             |      |      |      |
| 51   | 1   | CONNECTOR - HOUSING   | X-999046-02  |     |        |             |      |      |      |
| 52   | 1   | CONNECTOR - SOCKET  | X-999882-03  |     |        |             |      |      |      |
| 53   |     |   |  |     |        |             |      |      |      |
| 54   |     |   |  |     |        |             |      |      |      |
| 55   |     |   |  |     |        |             |      |      |      |

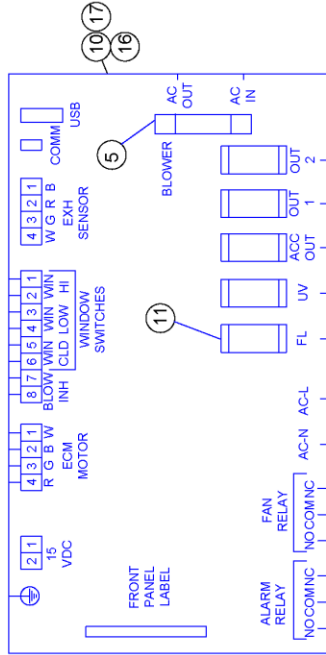

TITLE  
ELECTRICAL SCHEMATIC NU-240 115V 60 HZ  
(SERIES B)

|                |    |      |      |    |       |           |    |   |
|----------------|----|------|------|----|-------|-----------|----|---|
| DFTR           | TH | DATE | CHKD | BK | SHEET | 2         | OF | 2 |
| DRAWING NUMBER |    |      |      |    |       | CD-000098 |    |   |
|                |    |      |      |    |       | E         |    |   |