NU-543E CLASS II TYPE A2 BIOSAFETY CABINET



NJARE ZA




FlowGard<sup>™</sup> Electronic Control System

The LabGard<sup>®</sup> ES (Energy Saver) AIR Model NU-543E Class II, Type A2 Biosafety Cabinet (BSC) is a bench/ table top model with available base stand options. The NU-543E offers economical operation through the use of our single ultra-high efficiency DC ECM blower to reduce energy consumption and maximize filter load capacity. The NU-543E also offers the greatest degree of versatility of available options and accessories as well as a 254mm access opening. As a Class II Type A2, the NU-543E can be recirculated back into the room or configured with a canopy transition to an external exhaust blower. The The FlowGard<sup>™</sup> central control system utilizes a single airflow sensor with dual-thermistor probes to monitor and display the airflow speed in real-time.

DESAIR



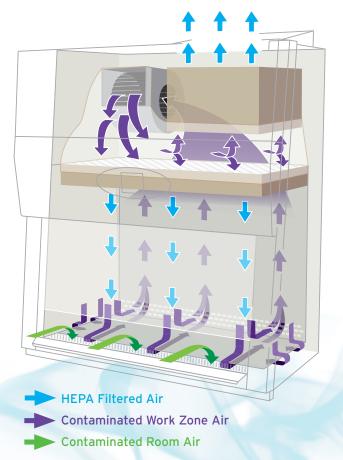

NuAire, Inc. | 2100 Fernbrook Lane | Plymouth, MN 55447 | U.S.A. | www.nuaire.com

## Popular Features & Accessories



	NU-543E Features & Specifications
Class, Type	Class II, Type A2
Cabinet Style	Bench Top w/ Optional Base Stand
Nominal Widths	0.9, 1.2, 1.5, 1.8 m
Workspace Dimensions (W x D x H)	(873 / 1178 / 1483 / 1788) x 654 x 724 mm
Exterior Dimensions (W x D x H)	(1057 / 1362 / 1669 / 1972) x 817 x 1546 mm
Net Weight	186 - 290 kg
Certifications	NSF / ANSI 49, EN12469, CE
Access Opening	254 mm
Electrical Requirements	230 VAC / 50-60 Hz
Filter Load Capacity	250%
Airflow	0.53 m/s Dynamic Inflow, 0.30 m/s True Laminar Downflow
	ISO Class 5 Air Protection
Supply Filter	76 mm H14 99.995% Efficient on MPPS
Exhaust Filter	292 mm H14 99.995% Efficient on MPPS
	Double Exhaust Filter
Control System	FlowGard™ with Audible/Visual Alarm
Airflow Sensor	IntelliFlow™ Airflow Sensor
Lighting	LED
	Ultraviolet (UV) Light Option
Nitecare™ Compatible	Yes
Password Protection	Yes
Decontamination Mode	Decon101™
Noise Level	Up to 63 dba*
Construction	Removable Stainless Steel Coved Work Tray with Prop-Up Stand
	2 Back-Wall Electrical Outlets
	3 Service Couplings and 1 Service Valve (Right/Left Sidewalls)
	Single Cord Pass-Through Port (Right Sidewall; NSF-Approved)
	Coved Interior Corners for Easy Cleaning
	Monolithic 16/18 gauge Stainless Steel Welded Pressure Tight Construction
	Recessed & Removable Front Intake Grill
	Stainless Steel Spill Trough w/ Drain Valve
	Attenumont™ Vibration Control Support System
	Permanent Positive Pressure Plenum w/ Quick Release Supply Filter Removal
	Frameless Easy-to-Slide Viewing Window w/ Optional Stainless Steel Edge Protector
	Heavy Duty Plastic Armrest w/ Stainless Steel Option
	IV Bar w/ 3 Height Locations, 6 SST Hooks
	Glass Sidewalls
	Backwall Cutout for PC Viewing Monitor
	3/8" Soft Copper, 3/8" Black Pipe (Gas Only), or 1/2" Rigid Copper Plumbing
	Additional Service Valves, Couplings, Cord Pass-Through Ports, and/or Outlets
Base Stand Accessories	Telescoping Adjustable Base Stand w/ Leg Levelers or Casters
	Motorized Adjustable Base Stand w/ Leg Levelers of Casters
Accessories	
	Variable Flow or Fixed Exhaust Canopy Transition
	Adjustable Ergonomic Footrest
	Lay-In Charcoal/Purifil Transition Options
	203 or 305 mm Base Stand Storage Shelves (Telescoping Only)
	Ergonomic Stainless Steel Turntable
	Silicone Cushioned Movable Elbow Rests
Standard Featu	re Optional Feature

\*Margin of error: ±2 dbA, measurement performed per ISO 11201:2010 in normal running mode

To Learn More or Request a Quote:

Visit: www.nuaire.com or Call: 1-763-553-1270 Toll Free: 1-800-328-3352



\*Contact your local NuAire distributor for international warranty details.













20-1656-F-I-EN-1-0923