NU-543S CLASS II TYPE A2 BIOSAFETY CABINET



7

NUAR




FlowGard[™] Electronic Control System

The LabGard[®] ES (Energy Saver) AIR Model NU-543S Class II, Type A2 Biosafety Cabinet (BSC) is a bench/ table top model with available base stand options. The NU-543S offers economical operation through the use of our single ultra-high efficiency DC ECM blower to reduce energy consumption and maximize filter load capacity. The NU-543S also offers the greatest degree of versatility of available options, accessories, and an access opening of 254 mm. As a class II Type A2, the NU-543S can be recirculated back into the room or configured with a canopy transition to an external exhaust blower. The FlowGard[™] central control system utilizes a single airflow sensor with dual-thermistor probe to monitor, display, and control airflow velocity in real-time.

FSAIR



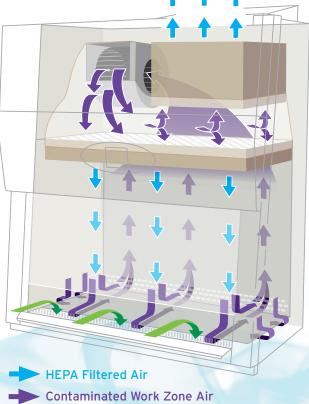

Contaminated Room Air

NuAire, Inc. | 2100 Fernbrook Lane | Plymouth, MN 55447 | U.S.A. | www.nuaire.com

Popular Features & Accessories



NU-543S Features & Specifications	
Class, Type	Class II, Type A2
Cabinet Style	Bench Top w/ Optional Base Stand
Nominal Widths	0.9, 1.2, 1.5, 1.8 m
Workspace Dimensions (W x D x H)	(873 / 1178 / 1483 / 1788) x 654 x 724 mm
Exterior Dimensions (W x D x H)	(1057 / 1362 / 1669 / 1972) x 817 x 1546 mm
Net Weight	186 - 290 kg
Certifications	NSF / ANSI 49, CE, UL listed
Access Opening	254 mm
Electrical Requirements	115V / 60Hz, 15-20 Amp
Filter Load Capacity	250%
Airflow	0.53 m/s Dynamic Inflow, 0.30 m/s True Laminar Downflow
	ISO Class 5 Air Protection
Supply Filter	76 mm 99.995% Efficient @ MPPS
Exhaust Filter	305 mm 99.995% Efficient @ MPPS
Control System	FlowGard™ with Audible/Visual Alarm
Airflow Sensor	IntelliFlow™ Airflow Sensor
Lighting	LED
	Ultraviolet (UV) Light Option
Nitecare™ Compatible	Yes
Password Protection	Yes
Decontamination Mode	Decon101™
Noise Level	Up to 63 dba*
Construction	Removable Stainless Steel Coved Work Tray with Prop-Up Stand
	Qty 2 Backwall Electrical Outlets
	3 Service Couplings and 1 Service Valve (Right/Left Sidewalls)
	Single Cord Pass-Through Port (Right Sidewall; NSF-Approved)
	Coved Interior Corners for Easy Cleaning
	Seam-Welded Stainless Steel Side and Back Wall Construction
	Recessed & Removable Front Intake Grill
	Stainless Steel Spill Trough w/ Drain Valve
	Attenumont [™] Vibration Control Support System
	Permanent Positive Pressure Plenum w/ Quick Release Supply Filter Removal
	Frameless Easy-to-Slide Viewing Window
	Heavy Duty Plastic Armrest
	IV Bar w/ 6 SST Hooks
	Additional Service Valves
Base Stand Accessories	Telescoping Adjustable Base Stand w/ Leg Levelers or Casters
Accessories	Fixed Exhaust Canopy Transition
Standard Feature Optional Feature	

*Margin of error: ±2 dbA, measurement performed per ISO 11201:2010 in normal running mode

To Learn More or Request a Quote:

Visit: www.nuaire.com or Call: 1-763-553-1270 Toll Free: 1-800-328-3352



*Contact your local NuAire distributor for international warranty details.









TIV NORD